

QUICK RESULTS, EASY TO USE, ACCURATE AND RELIABLE

On-site Blood Chemistry Analyzer

1. Add sample











EASY TO USE

Fully Automated System - no special operating skills required

- 8.4 inch capacitive touch screen, more friendly interface.
- 0.1cc whole blood, serum or plasma
- Barcoded prefabricated calibration information
- Built in printer, results printed directly
- LIS comaptible
- Results could be printed through an office printer or through data management platform (MMDMP)

ACCURATE AND RELIABLE

Advanced technology ensures precise results

- Microfluidics discs with pre-installed reagent pellets to
- ensure accurate analysis of blood samples and reagents
- Stable measurement optics include a troboscopic xenon
- lamp, a wavelength selection system, and a multiple-wavelength detector
- Integrated quality control software monitors the entire process in real time.

MUREX Diagnostic Products Specialists -










Celecare cM5 Reagent Discs

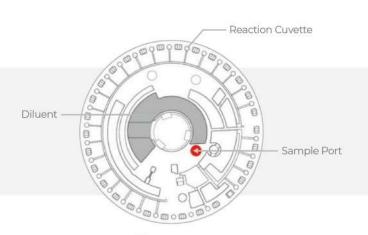

A 6	and the same	on Parish	Rapide C. Ashell .d.	Cardier Ath
See Schille St. Chica Schille See	Action of the Parties	the day of hitself and a	du Hadilled & Cherry Childed	and the state of t
Care of the State of the State of State	South of the South	Georgia Parille Parille	diland of the state of the stat	A Separate of the separate of

TP ALB GLO* A/G* ALT AST GGT ALP TBIL DBIL IBIL* UA CRE UREA CK CK-MB LDH a-HBDH AMY GLU GSP HCY TG CHOL HDL-C LDL-C* Na* CI Ca2*

00 3	CH. 74	50. A.	24 A	4. 14	Che A	0 4	Q 3	Co 7	2, 2	bo 12	Co. 7
TP			TP						TP		TP
ALB		ALB	ALB						ALB		ALB
GLO"			GLO"						GLO*		GLO*
A/G*			A/G"						A/G*		A/G*
ALT			ALT						ALT		ALT
AST	AST		AST	AST					AST		AST
			GGT						GGT		GGT
			ALP								ALP
TBIL			TBIL						TBIL		TBIL
DBIL			DBIL								DBIL
IBIL*			IBIL*								IBIL*
UA	UA	UA									UA
CRE	CRE	CRE						CRE	CRE		CRE
		UREA						UREA	UREA		UREA
	CK			CK							
	CK-MB			CK-MB							
	LDH			LDH							
	a-HBDH			а-НВОН							
	AMY							AMY			
GLU	GLU					GLU	GLU	GLU	GLU		GLU
						GSP					
							HCY				
TG						TG	TG				TG
CHOL						CHOL	CHOL				CHOL
HDL-C						HDL-C	HDL-C				HDL-C
LDL-C*						LDL-C*	LDL-C*				LDL-C*
	K.				K.			K'			
	Na*				Na*			Na*			
	CI				CI			CI			
		Ca2*			Ca2*						
		P			Р						
					Mg ²⁺						
	CO2	CO2			CO			CO,			
										NH,	

* Calculated test value

MUREX Diagnostic Products Specialists -



P Mg²⁺ CO₂ NH₃



