



## **BC-30s**

Auto Hematology Analyzer

Minimum Size,  
Maximum Capability



**mindray**  
healthcare within reach

# BC-30s


## Auto Hematology Analyzer

### New 3-part

With decades' dedication and cumulative experience in establishing hematology analyzer platform, BC-30s is newly designed for diagnosis labs wishing to have better usability and higher efficiency.

BC-30s, the minimum size hematology analyzer with extraordinary exquisite appearance Mindray have ever been designed, may offer you advanced user friendly operation, detailed flag information, which is the most compact footprint and the best power efficiency 3-part hematology analyzer with maximum capability.

Parameter	Result	Unit
WBC	L 2.4	10 <sup>9</sup> /L
Lym#	L 0.7	10 <sup>9</sup> /L
Mid#	R 0.3	10 <sup>9</sup> /L
Gran#	RL 1.4	10 <sup>9</sup> /L
Lym%	0.305	
Mid%	R 0.130	
Gran%	R 0.565	

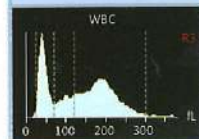

Leucopenia  
Lymphocyte  
decreased  
Pancytopenia

## Exclusive Feature



The screenshot displays the BC-30s interface with a detailed flag message for WBC. The message is highlighted in a red box and reads:

- WBC Message
- Lymphocyte increased
- RBC Message
- RBC Distribution Abnormal
- Anemia
- HGB Abn./Interfere?
- PLT Message

The interface also shows a table of hematology parameters:

Parameter	Result	Unit	Parameter	Result	Unit	Parameter	Result	Unit
WBC	RH 10.9	10 <sup>9</sup> /L	RBC	L 1.10	10 <sup>12</sup> /L	PLT	L 74	10 <sup>9</sup> /L
Lym#	RH 7.9	10 <sup>9</sup> /L	HGB	RL 51	g/L	MPV	H 13.0	fL
Mid#	R 0.8	10 <sup>9</sup> /L	HCT	L 0.110		PDW	H 17.9	
Gran#	R 2.3	10 <sup>9</sup> /L	MCV	H 100.3	fL	PCT	L 0.96	mL/L
Lym%	RH 0.723		MCH	RH 46.8	pg	P-LCC	37	10 <sup>9</sup> /L
Mid%	R 0.069		MCHC	RH 466	g/L	P-LCR	H 0.495	
Gran%	RL 0.208		RDW-CV	H 0.182				
			RDW-SD	H 68.1	fL			

Detailed flag message, as exclusive feature, assists doctors for pre-diagnosics. It can provide flag information which are similar with 5-part analyzer including WBC flag, RBC flag and PLT flag

Different flag information are provided according to parameters results together with histograms.



**Sample 1** : BC-30s can provide flags named "Leucopenia", "Lymphocyte decreased", "Pancytopenia" which represent white blood cell decreased, the low number of Lymphocyte and decreased of leukocyte, erythrocyte and plate count. Also the previous flag "R3" is still shown. Two kinds of flag messages are both supported to ensure clinicians have better understanding of sample results.

**Sample 2** : Flag "Anemia" is pop up. It means the sample has the possibility of anemia.

**Sample 3** : Flag "Thrombopenia" is displayed which stands for Plates decreased.


Parameter	Result	Unit
RBC	L 1.86	10 <sup>12</sup> /L
HGB	L 67	g/L
HCT	L 0.193	
MCV	H 103.5	fL
MCH	H 35.8	pg
MCHC	346	g/L
RDW-CV	0.157	
RDW-SD	H 60.7	fL

Anemia

Parameter	Result	Unit
PLT	L 38	10 <sup>9</sup> /L
MPV	H 12.6	fL
PDW	H 18.5	
PCT	L 0.48	mL/L
P-LCC	L 17	10 <sup>9</sup> /L
P-LCR	H 0.455	



Thrombopenia

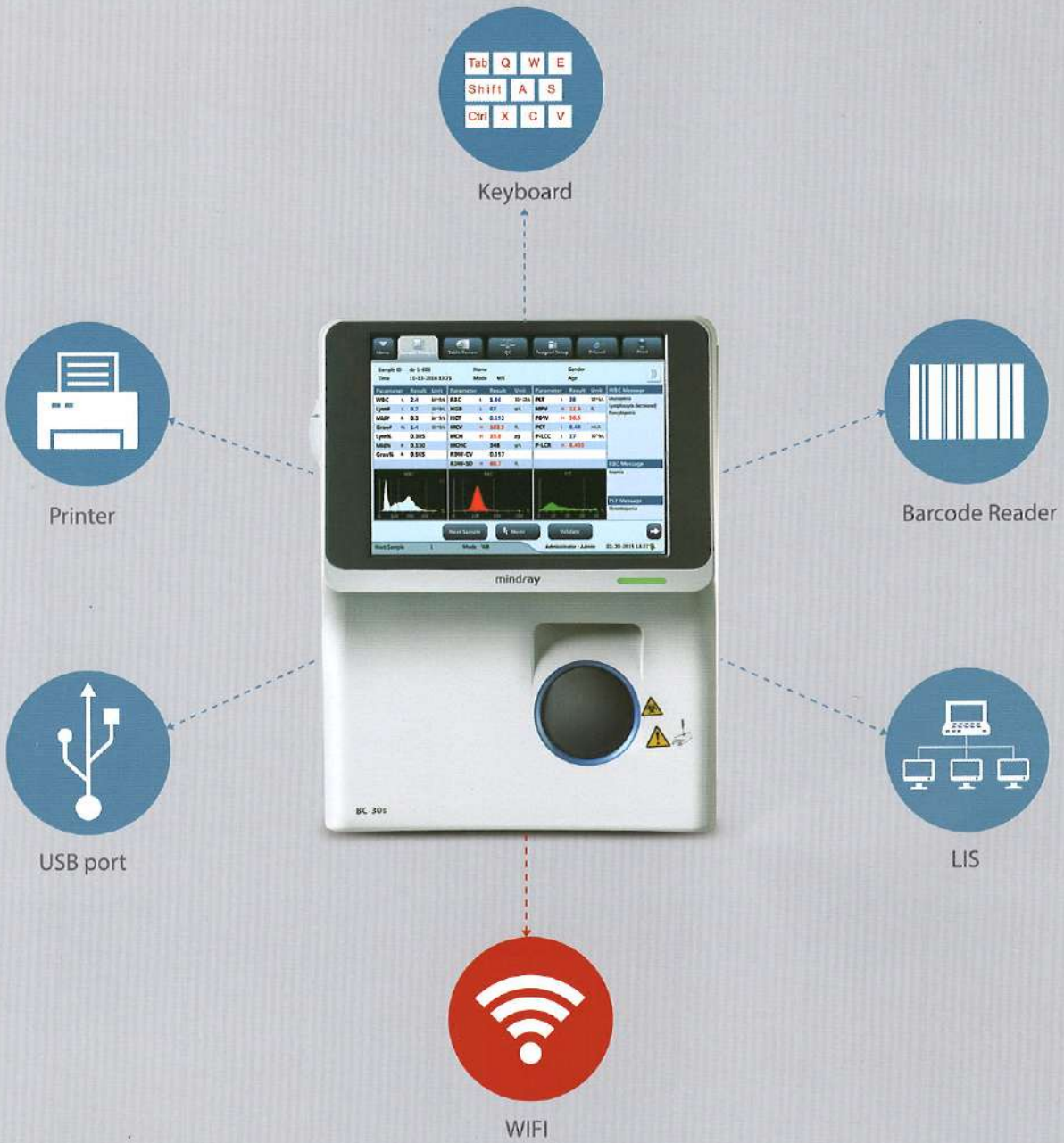

WIFI capability provides you full solution for data communication together with bi-directional LIS, USB port and LAN port, barcode reader, printer and keyboard.



## Better Usability



### Slim size spares you space

Minimum size with small footprint of 17 inch laptop as well as Internally-stored-lyse design save more space for small labs to ensure them have higher-usage of room.

10.4 inch TFT touch screen together with powerful software offer you more user-friendly interface and pleasant operation experience.



### Smart design offers you lower reagent consumption

Less reagents is needed and lower reagent consumption is effected owing to optimized liquid system, EZ cleanser and rinse deleted.

In order to match the requirements of low sample volume target customers, smaller package of reagents are supported.

## Enhanced Quality

### Swift speed provides you good specification

Higher throughput at 70 tests per hour, it will save clinicians more time.

### Superior features ensure you higher performance

Micro sample volume at 9.0uL for whole blood mode along with capillary whole blood samples supported are perfect for pediatric samples.



# BC-30s

## Auto Hematology Analyzer

### Technical Specifications

#### Principles

Impedance method for WBC, RBC and PLT counting  
Cyanide free reagent for hemoglobin test

#### Performance

Parameter	Linearity Range	Precision (CV %)	Carryover
WBC( $10^9/L$ )	0-200	$\leq 3.5\%$ (4.0-6.9) $\leq 2.0\%$ (7.0-15.0)	$\leq 0.5\%$
RBC( $10^{12}/L$ )	0-8.00	$\leq 1.5\%$ (3.5-6.5)	$\leq 0.5\%$
HGB(g/L)	0-280	$\leq 1.5\%$ (100-180)	$\leq 0.5\%$
MCV(fL)		$\leq 1.0\%$ (70-110)	
PLT( $10^9/L$ )	0-4000	$\leq 5.0\%$ (100-149) $\leq 4.0\%$ (150-500)	$\leq 1.0\%$

#### Parameters

21 parameters: WBC, Lymph#, Mid#, Gran#, Lymph%, Mid%, Gran%, RBC, HGB, HCT, MCV, MCH, MCHC, RDW-CV, RDW-SD, PLT, MPV, PDW, PCT, P-LCR, P-LCC  
3 histograms for WBC, RBC and PLT

#### Reagent

M-30D DILUENT  
M-30CFL LYSE  
PROBE CLEANSER

#### Sample Volume

Prediluted mode 20 $\mu$ L  
Whole blood mode 9 $\mu$ L

#### Throughput

70 samples per hour

#### Display

10.4 inch TFT Touch Screen

#### Data Storage Capacity

Up to 500,000 results including numeric and graphical information

#### Communication

LAN Port supports HL7 protocol  
Support bi-directional LIS

#### Interface

4 USB port (for external printer, software upgrade, barcode reader, WIFI adapter, keyboard and mouse), LAN port (1)

#### Printout

Thermal recorder, 50 mm wide paper, various printouts formats  
External printer optional

#### Operating Environment

Temperature: 10°C~40°C  
Humidity: 10%~90%  
Air pressure: 70kPa~106kPa

#### Power Requirement

100V-240V  
 $\leq 300VA$   
50Hz/60Hz

#### Dimension and Weight

Dimension: Depth(410 mm) x width(300 mm) x height(400 mm)  
Weight:  $\leq 20Kg$

Mindray Building, Keji 12th Road South,  
High-tech Industrial Park, Nanshan, Shenzhen 518057, P.R. China  
Tel: +86 755 8188 8998 Fax: +86 755 26582680  
E-mail: intl-market@mindray.com www.mindray.com

**Mindray is listed on the NYSE under the symbol "MR"**

**mindray** are registered trademarks or trademarks owned by Shenzhen Mindray Bio-medical Electronics Co., Ltd.  
© 2015 Shenzhen Mindray Bio-Medical Electronics Co., Ltd. All rights reserved. Specifications subject to changes without prior notice.  
P/N: ENG-BC-30s-210285x6-20141229

**mindray**



# **BC-5390**

## **Auto Hematology Analyzer**



**mindray**  
healthcare within reach

# BC-5390


## Auto Hematology Analyzer

The new BC-5390 Auto Hematology Analyzer provides rapid and reliable test from just 33uL of blood. Utilizing three mainstream technologies: laser scatter, flow cytometry and chemical dye, BC-5390 can provide accurate differential readout. In order to save time and increase walk-away automation, an autoloader is equipped to hold 40 tubes one time and achieve 60 samples/hr throughput. With the popular windows based software, you can easily perform the routine tests, manage patient results, set up auto-cleaning and connect with LIS server. BC-5390 is your ideal choice to streamline daily workflow.




- Semi-conductor laser combined with chemical dye method advanced flow cytometry
- 21 parameters with complete 5-part differentiation(CBC+DIFF) on white blood cells
- 40-tubes autoloader with random access
- Closed tube for STAT samples
- Capability to flag abnormal samples
- Only 33  $\mu$ l sample volume for CBC + DIFF results
- Up to 60 samples processed per hour
- Supports whole blood mode for capillary sample
- Large storage capacity: up to 100,000 samples
- Support bi-directional LIS connection
- Customization on reference range, auto-cleaning schedule and report format





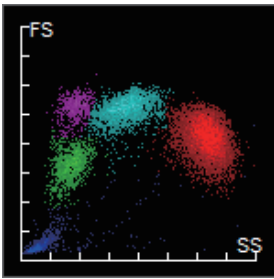

Eosinophilia sample

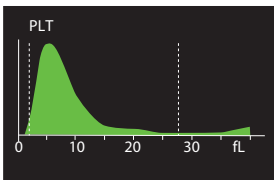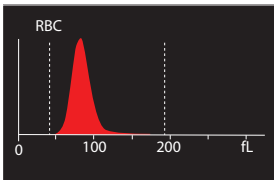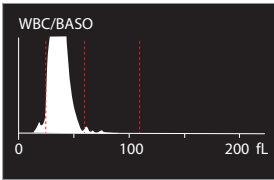

**White blood cell differentiation**

The semi-conductor laser system collects forward scatter and side scatter information of white blood cells in a flow cell and differentiates lymphocytes, monocytes, neutrophils and eosinophils according to the cell size and granule complexity. Specific chemical dye treatment to eosinophils can separate from neutrophils in DIFF scattergram even on Eosinophilia sample . In addition, basophils are measured in an impedance channel upon lysing action on the RBC and other WBC groups.


**Advanced impedance measurement for RBC and PLT**

Mindray patented digital sheath flow and weak signal detection techniques are applied to ensure accurate counting and sizing of RBC and PLT in the impedance channel. Dynamic discriminator adjustment can improve the separation of RBC and PLT groups when large platelet and micro RBC are present. Benefiting from these improvements, MCV, RDW and PLT performance are greatly enhanced.



**M-5 reagents and BC-5D controls**

Only 4 routine M-5 reagents are utilized in counting cycle featuring economical consumption and 2 years shelf life. To help precision monitoring, three levels of BC-5D controls are offered in a ready-to-use kit and the assay value table can be automatically imported through USB memory.

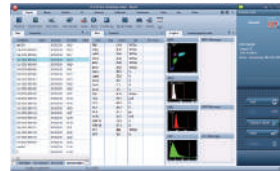

**Windows based software**

The analyzer's windows-based software is simple to use and plays a powerful information hub to store 100000 patient results. You can set-up password access, reference range, auto-cleaning schedule and so on. Also, the built-in report format tool can help to customize the final report type to include microscopic counting, ESR, blood type and diagnostic remarks.




**QC monitoring and patient archive**

60 QC files are designed to store L-J QC results. 300 data points can be recorded in single file. 2 common QC programs are enabled for full quality assurance purposes. Patient data is archived and can be searched and presented in trend curve for case follow-up.



Venous whole blood sample    Pre-diluted sample    Capillary whole blood sample

Different holder adapters supports three types of sample including Venous whole blood sample, Predilute sample and Capillary whole blood sample. Capillary whole blood sample can be well tested directly which is more convenient for users and makes the analyzers an ideal choice to work with pediatric and geriatric samples.

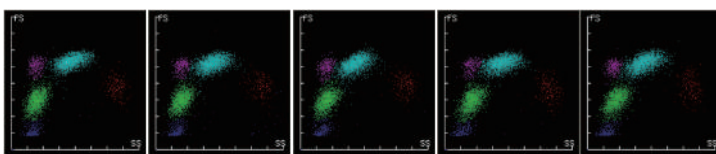

1 hour    4 hours    8 hours    12 hours    24 hours

Anti-aging samples with more than 12 hours transportation to central labs which is a quite frequent and normal situation for small to medium labs in countryside. Hence, diagnosing anti-aging samples accurately becomes an significant issue for developing countries or cities. With enhanced reagent system and upgraded reagent formula, BC-5390 is an ideal solution that can better differentiate white blood cells on anti-aging samples under room temperature.

# BC-5390

## Auto Hematology Analyzer

### Technical Specifications:

#### Principles

Impedance method for WBC/BAS, RBC and PLT counting;  
Cyanide free reagent for hemoglobin test;  
Flow Cytometry (FCM) + Laser scatter + Chemical dye method for  
WBC differential analysis

#### Parameters

21: WBC, LYM(#,%), NEU(#,%), MON (#,%), EOS (#,%), BAS (#,%), RBC, HGB,  
HCT, MCV, MCH, MCHC, RDW-CV, RDW-SD, PLT, MPV  
3 histograms and 1 scattergram

#### Analysis Mode

Autoloader Whole Blood Mode  
Closed Whole Blood Mode  
Closed Predilute Mode

#### Throughput

Up to 60 samples per hour with autoloader  
Up to 51 samples per hour with closed mode  
Up to 53 samples per hour with predilute mode

#### Sample Volume

Whole Blood CBC+DIFF 33  $\mu$ L  
Whole Blood CBC 24  $\mu$ L  
Prediluted: 20  $\mu$ L

#### Test Panel

CBC CBC+DIFF

#### Performance

	Carryover	Precision	Linearity
WBC ( $\times 10^3/\mu$ L)	$\leq 1.0\%$	$\leq 0.15$ (SD) or 3.0% (CV%)	0.3-200
RBC ( $\times 10^6/\mu$ L)	$\leq 1.0\%$	$\leq 1.5\%$	0.2-8.0
HGB (g/dL)	$\leq 1.0\%$	$\leq 1.5\%$	0.5-25
HCT (%)	$\leq 1.0\%$	$\leq 2.0\%$	2 -75
PLT ( $\times 10^3/\mu$ L)	$\leq 1.0\%$	$\leq 7.5$ (SD) or 5.0% (CV%)	5-2000

#### Data Storage Capacity

Up to 100,000 results including numeric and  
graphical information

#### Communication

LAN Port supports HL7 protocol

#### Operating Environment

Temperature: 15  $^{\circ}$ C~30 $^{\circ}$ C  
Humidity: 30~85%  
Air Pressure: 70~106 kPa

#### Power Requirement

A.C.100-240V  $\leq$  300VA 50/60Hz

#### Dimension and Weight

Dimension (WxDxH, inches): 22.4 x 23.2 x 20.6  
Weight:  $\leq$  143 pounds

# BC-760 & BC-780

Auto Hematology Analyzer with ESR

## Key Specifications

### Principles

WBC (IMG/Neu/Mon/Lym/Eos/Bas), NRBC/RET\*, PLT-H/PLT-O\*/IPF:  
SF Cube ^ Cell Analysis Technology  
^S: Scatter; F: Fluorescence; Cube: 3D analysis

### RBC, PLT

Focusing Flow-DC Impedance Method

### HGB

Colorimetric method

### ESR

Photometric method

**Number of measuring parameters (whole blood): 109**

**Number of reportable parameters: 41**

WBC Bas# Bas% Neu# Neu% Eos# Eos% Lym# Lym% Mon#  
Mon% IMG# IMG% RET#\* RET#\* RHE\* IRF\* LFR\* MFR\* HFR\*  
RBC HGB MCV MCH MCHC RDW-CV RDW-SD HCT NRBC#  
NRBC% PLT PLT-I PLT-H PLT-O\* MPV PDW PCT P-LCR P-LCC  
IPF ESR

**Number of research parameters: 68\***

**Number of measuring parameters (body fluid): 18**

**Number of reportable parameters: 7**

WBC-BF TC-BF# MN# MN% PMN# PMN% RBC-BF

**Number of research parameters: 11**

### Sample volume

CD (whole blood): 25ul CD+ESR

(whole blood): 160ul Predilute: 20ul

### Data storage capacity

Up to 150,000 results including numeric and graphical information \*

### Throughput

CD 80t/h CDR 45t/h CD+ESR 40t/h

### Analysis Mode

Sample Type	Analysis Mode
Whole blood	CBC, CBC + DIFF, CBC + DIFF+RET*, CD + ESR, CDR + ESR*, CD/WBC-3X, CDR/PLT-5X*, and other modes
Predilute	CBC, CBC + DIFF, CDR*, and other modes
Body fluid	CBC + DIFF

## Physical Specifications

### Dimensions

840D x 655W x 600H mm

### Weight

≤70.6Kg

### Voltage

100V-240V~ (±10%)

### Frequency

50Hz/60Hz (±1Hz)

### Power input

600VA

### External output

LANx1, USB x 4 (Specifications: DC 5V; 500mA;

USB2.0 x 3; USB3.0 x 1)

## Normal Operating Environment

### Ambient temperature:

10°C ~ 35°C

### Relative humidity:

30% ~ 85%

### Atmospheric pressure:

70.0kPa ~ 106.0kPa^

^Note : Required altitude for normal operation:

-400m ~ +3000m

## Performance

Parameter	Linearity Range	Precision	Carryover
WBC	0-500×10 <sup>9</sup> /L	≤2.5% (≥4.51×10 <sup>9</sup> /L)	≤ 1.0%
RBC	0-8.60×10 <sup>12</sup> /L	≤1.5% (≥3.5×10 <sup>12</sup> /L)	≤ 1.0%
HGB	0-260g/L	≤1.0% (110-180g/L)	≤ 1.0%
HCT	0-75%	≤1.5% (30%-50%)	≤ 1.0%
PLT*	0-5000×10 <sup>9</sup> /L	≤ 1.5(SD) (≤20×10 <sup>9</sup> /L)^ ≤ 2.5% (≥100×10 <sup>9</sup> /L)^	≤ 1.0%
RET*	0-0.8×10 <sup>12</sup> /L	≤15% (RBC ≥ 3.00×10 <sup>12</sup> /L RET%: 1.00% ~ 4.00%)	≤ 1.0%
ESR		≤1.8(SD)(0~20mm/h)	≤ 1.0%

Note: Applicable only to CDR/PLT-O 5x and CR/PLT-O 5x models

Items marked with an asterisk (\*) apply only to BC-780

**mindray**

# BC-760 & BC-780

Auto Hematology Analyzer with ESR

## Above and Beyond



[www.mindray.com](http://www.mindray.com)

P/N:ENG-BC-760 & BC-780 ESR-210285X6P-20220115  
©2021 Shenzhen Mindray Bio-Medical Electronics Co., Ltd. All rights reserved.


**mindray**  
healthcare within reach

••• Above your expectations


► SF Cube fluorescent technology allows reliable counting and differentiation of abnormal samples

— More refined and reliable cell differentiation


3D fluorescent analysis technology allows reliable differentiation of immature and other abnormal cells, such as immature granulocytes (IMGs), reticulocytes (RETs\*), and immature platelet fraction (IPF).

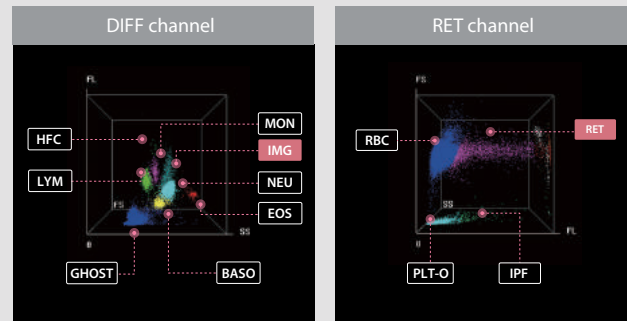

WBC-Part Differentiation

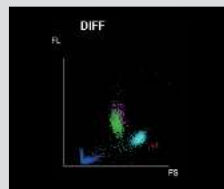


Conventional

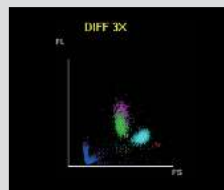


WBC 3x

PLT-O Measurement

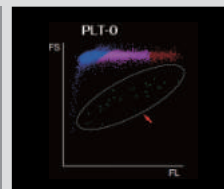

Conventional

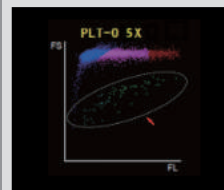

PLT-O 5x

— More reliable measurements for low-value samples

The BC-760 & BC-780 3D fluorescence analysis platform is designed with multiple counting WBC-3x and PLT-O 5x analysis modes to help ensure higher reliability for low-value WBC and PLT samples. In addition, the PLT de-aggregation function can reduce the cumbersome review work.

— More comprehensive alarm messages for abnormalities


The analyzer provides a detailed list of over 40 prompt messages, including WBC message, RBC message, and PLT message. This allows laboratory technicians to intuitively and quickly identify abnormal samples and proceed further with the samples in a timely manner. This in turn helps to avoid missed diagnosis of blood disease and false reports.

Parameter	Result	Unit	Parameter	Result	Unit
WBC	18.88	10 <sup>9</sup> /L	RBC	4.89	10 <sup>12</sup> /L
Neu#	15.90	10 <sup>9</sup> /L	HGB	L	100 g/L
Lym#	0.32	10 <sup>9</sup> /L	HCT	L	0.362
Mid#	0.43	10 <sup>9</sup> /L	MCV	L	72.4 fL
Eos#	0.01	10 <sup>9</sup> /L	MCH	L	21.7 pg
Bas#	0.02	10 <sup>9</sup> /L	MCHC	L	300 g/L
RDW	0.43	10 <sup>9</sup> /L	RDW-CV	H	2.298 %
Neu%	0.564	%	RDW-SD	H	77.7 fL
Lym%	0.018	%	RETs	0.0304	10 <sup>12</sup> /L
Mid%	0.028	%	RET%	0.61	%
Eos%	0.000	%	IPF	10.6	%
Bas%	0.001	%	LRF	99.4	%
MPV	0.028	fL	MFR	8.0	%
PLT	4.1	10 <sup>9</sup> /L	IFR	2.6	%
MPV	10.8	fL	RW	L	25.5 pg
PDW	14.8	%	NRBC	0.000	10 <sup>12</sup> /L
PCF	1.31	%	NRBC%	0.00	%
PLCR	39	10 <sup>9</sup> /L	ESR	5.37	mm/h
IPF	45.3	%			
IPF	5.7	%			

••• Beyond your expectations


► Limitations of traditional PLT counting

In the traditional impedance method, PLTs are subject to interferences that may lead to falsely high or falsely low results (as shown in the figure). Once an error report is generated, it will directly affect the judgment and decision-making of clinicians. The results reported at the clinical decision level are related to patient safety. Therefore, accurate PLT results are critical in clinical practice.

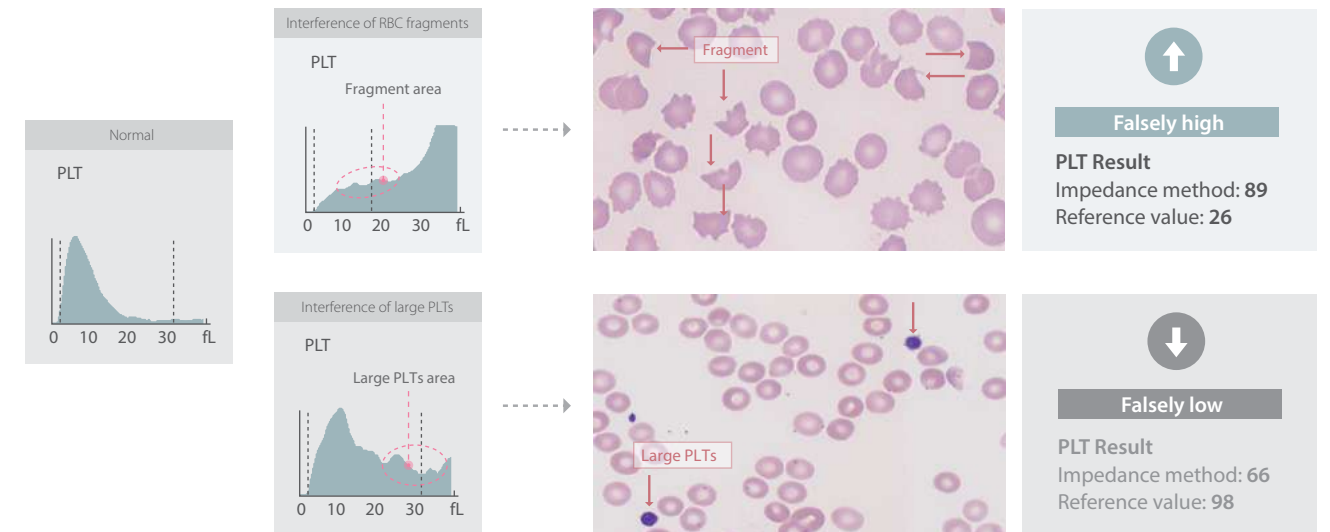

Optical PLT-H in every CD test

In order to solve the above problem, we have developed a brand new parameter PLT-H. It combines small PLTs from the conventional impedance method and large PLTs from the optical method. The solution can resist the interferences in conventional PLT detection without requiring extra reagents.

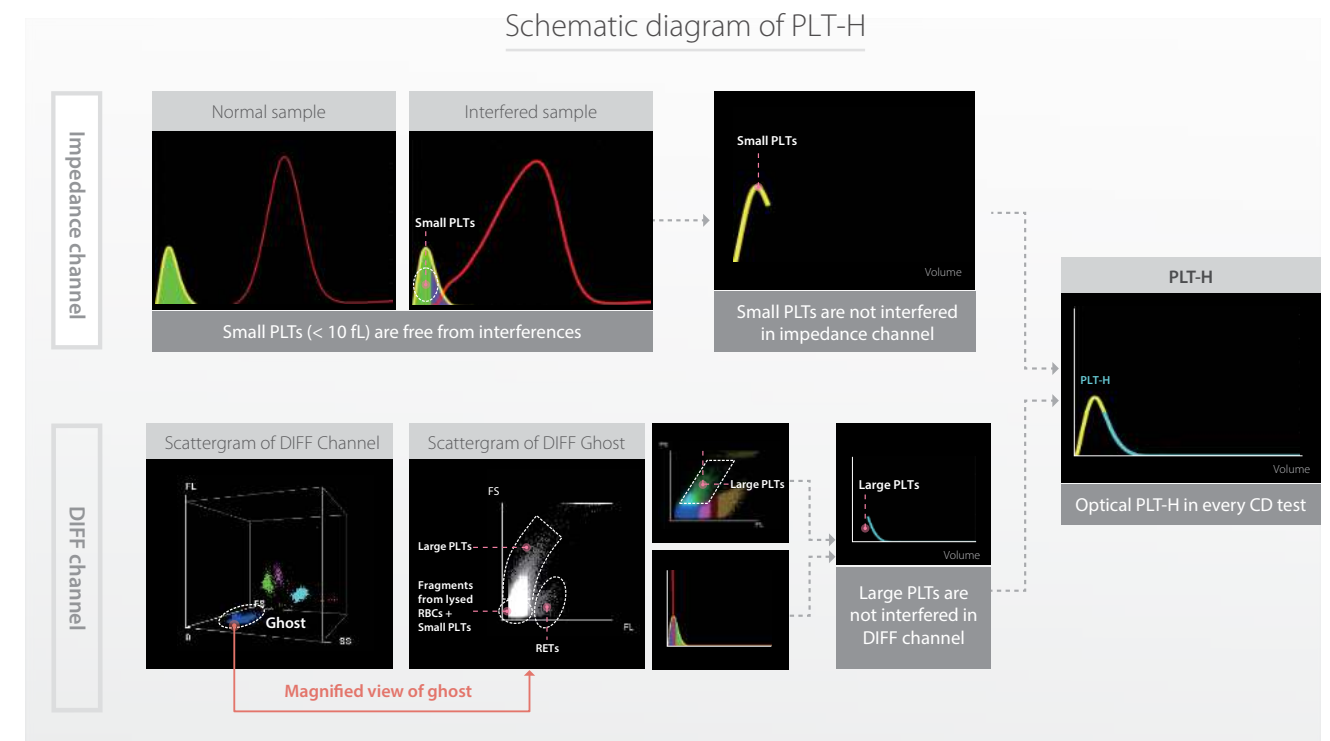



**RFID**  
Encryption key management



**labXpert**  
Coming as standard configuration  
Same software as BC-6 series



**Floating screen**  
Switching between different  
analysis modes with a single touch


▶ Excellent performance, high reliability, and ease of use



**Continuous autoloading**  
5 positions x 6 racks



**Cap piercing STAT**  
Supports STAT samples and capillary blood samples;  
minimizes bio-safety hazards

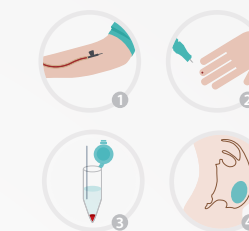

**Applicable to different sample types**  
Peripheral blood/Capillary  
blood Pre-diluted blood/Body  
fluid

▶ An all-in-one solution that goes above and beyond your expectations

## CD + ESR in one test provide reliable ESR results with greater ease

The BC-700 series integrates an automatic ESR module in a hematology analyzer. It can also generate both CBC & ESR results in one test within 1.5 min. In addition, it saves the costs that would otherwise be incurred for the purchase, maintenance, consumables, and storage space of a separate ESR analyzer. Compared with the traditional Westergren method, this method performs better in quality traceability, repeatability, speed, safety, and level of automation.

### Accurate

- Great correlation with the Westergren method
- Same QC and calibrator as in the BC-6000 series
- Combined examination helps to avoid the interferences of dehydration, polycythemia vera and anemia on ESR results



### Cost-effective

- The integrated instrument is capable of both CBC and ESR detection;
- Takes up the space of only one analyzer.



### Automatic

- Report CBC + ESR results together within 1.5 min;
- The measurement results are protected against the influence of subjective factors;
- Automation can reduce the biosafety hazards that may otherwise be introduced by a manual method.



### Excellent performance

Reliable results for abnormal and interfered samples

### SF Cube

- SF Cube 3D fluorescent platform
- IMG/NRBC reportable parameters
- Accurate measurement of low WBC/PLT count
- Anti-aging capability:  
24h (room temperature)  
48h (refrigerated transportation)

### Superior functionality


Integrated CBC + ESR detection provides a comprehensive tool for inflammation detection  
Optical PLT-H ensures reliable PLT results even with interference

- Integrated CBC + ESR analyzer
- Optical PLT-H in every CD test
- Integrated RET\* & body fluid analyzer

### Interactive excellence

Refined details for the optimal user experience

- iHelp videos
- Floating screen
- Ease of maintenance

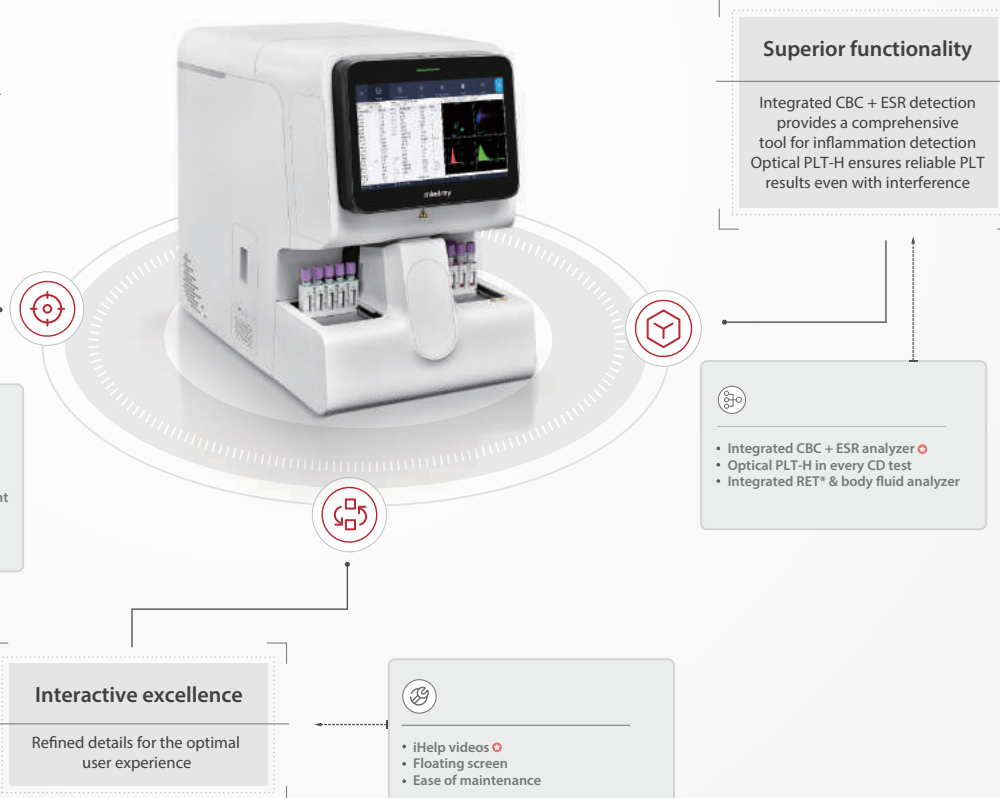

**Principles**

SF Cube\* method to count WBC, 6-Part diff and NRBC  
 DC impedance method for RBC and PLT  
 Cyanide free reagent for hemoglobin test  
 \*S: Scatter; F: Fluorescence; Cube: 3D analysis

**Parameters**

29 Reportable parameters (whole blood): WBC, Lym%, Mon%, Neu%, Bas%, Eos%, IMG%, Lym#, Mon#, Neu#, Eos#, Bas#, IMG#;  
 RBC, HGB, HCT, MCV, MCH, MCHC, RDW-CV, RDW-SD, NRBC#, NRBC%;  
 PLT, MPV, PDW, PCT, P-LCR, P-LCC  
 20 Research parameters (whole blood): HFC#, HFC%, WBC-D, TNC-D, IME%, IME#, H-NR%, L-NR%, NLR, PLR, WBC-N, TNC-N, InR#, InR%, Micro#, Micro%, Macro#, Macro%, PDW-SD, PLT-I  
 7 Reportable parameters (body fluid): WBC-BF, TC-BF#, MN#, MN%, PMN#, PMN%, RBC-BF  
 11 Research parameters (body fluid): Eos-BF#, Eos-BF%, Neu-BF#, Neu-BF%, HF-BF#, HF-BF%, RBC-BF, LY-BF#, LY-BF%, MO-BF#, MO-BF%

2 Histograms for RBC and PLT  
 2 Three-dimension scatter grams: DIFF, WNB  
 2 Two-dimension scatter grams: DIFF, WNB

**Mode**

CBC, CBC+DIFF

**Data storage capacity**

Up to 10,0000 results including numeric and graphical information

**Operating environment**

Temperature: 15 C ~32 C  
 Humidity: 30%~85%

**Performance**

Parameter	Linearity Range	Precision	Carryover
WBC	0-500×10 <sup>9</sup> /L	≤2.5% (≥4×10 <sup>9</sup> /L)	≤1.0%
RBC	0-8.60×10 <sup>12</sup> /L	≤1.5% (≥3.5×10 <sup>12</sup> /L)	≤1.0%
HGB	0-260g/L	≤1.0% (110-180g/L)	≤1.0%
HCT	0-75%	≤1.5% (30%-50%)	≤1.0%
PLT	0-5000×10 <sup>9</sup> /L	≤4.0% (≥100×10 <sup>9</sup> /L)	≤1.0%

**Sample volume**


Whole blood (Autoloader, Closed Tube)	80uL
Capillary blood (Closed Tube)	35uL
Predilute (Closed Tube)	20uL
Body fluid (Closed Tube)	85uL

**Throughput**

Up to 110 samples per hour (CBC+DIFF)  
 Up to 40 samples per hour (Body fluid)

**Loading capacity**

Up to 50 sample tubes

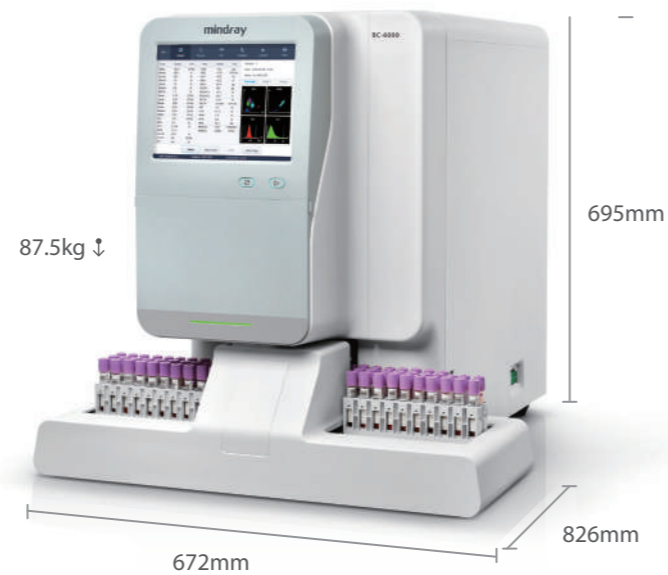

**BC-6000**

Auto Hematology Analyzer


**High Performance for ALL**






At Mindray, we seek to understand the needs of every customer, and deliver tailor-made solutions. Before designing any product, we listen to the **Voice Of Customers** and bear in mind the challenges they face. In today's laboratories, lab managers are looking for an analyzer with greater clinical values, such as higher flagging efficiency to reduce the ratio of microscopic examination, NRBC/body fluid results generated in a small-footprint system, among others, all within limited budget.

Mindray BC-6000 fulfills all these requirements and exceeds your expectations with so much more in terms of **Value**, **Operation** and **Cost**.

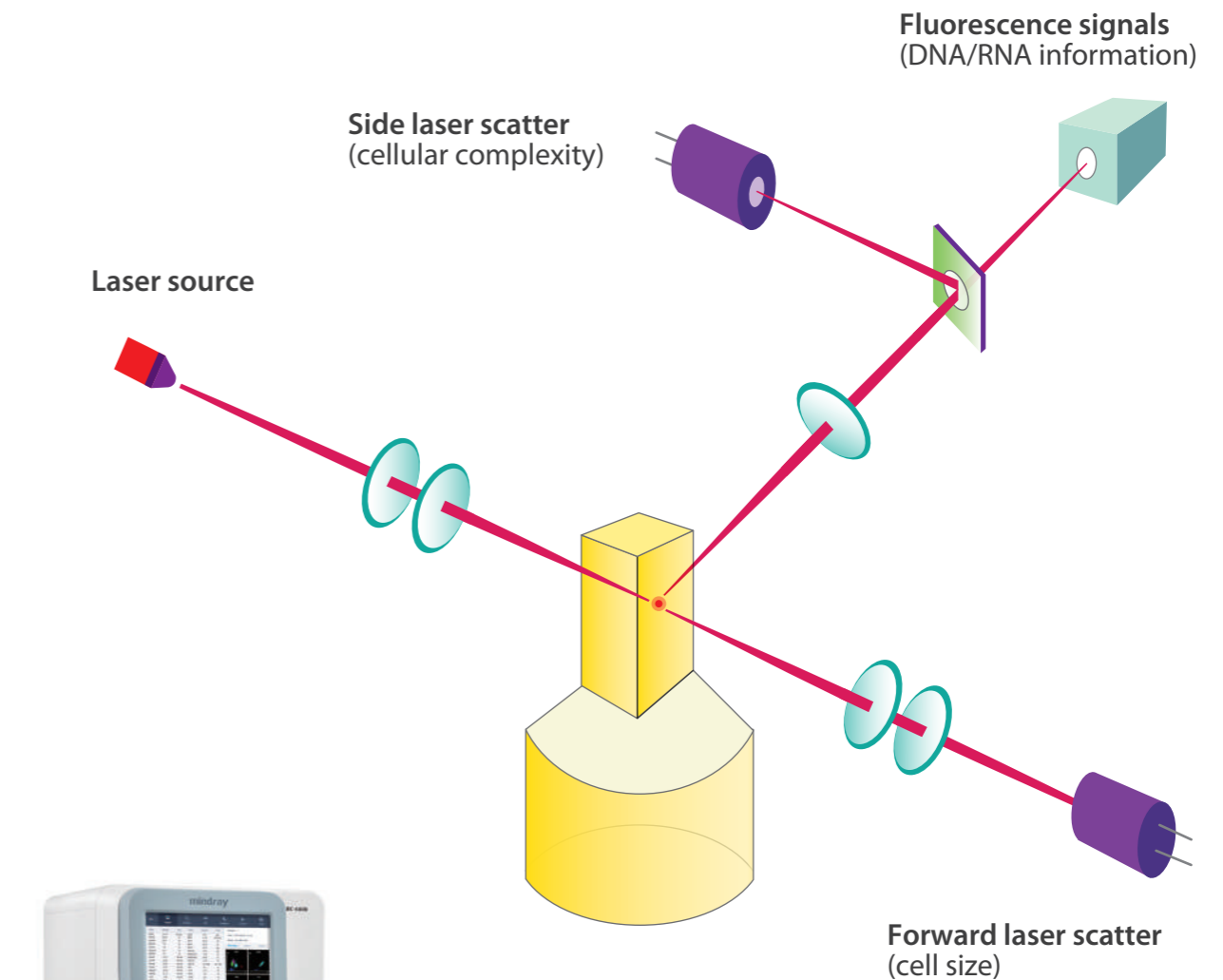

With the newly designed optics and reagent systems, the SF Cube technology can help doctors to better differentiate the clusters of cells, which is the key to revealing more abnormal cells.


# DIFF Channel

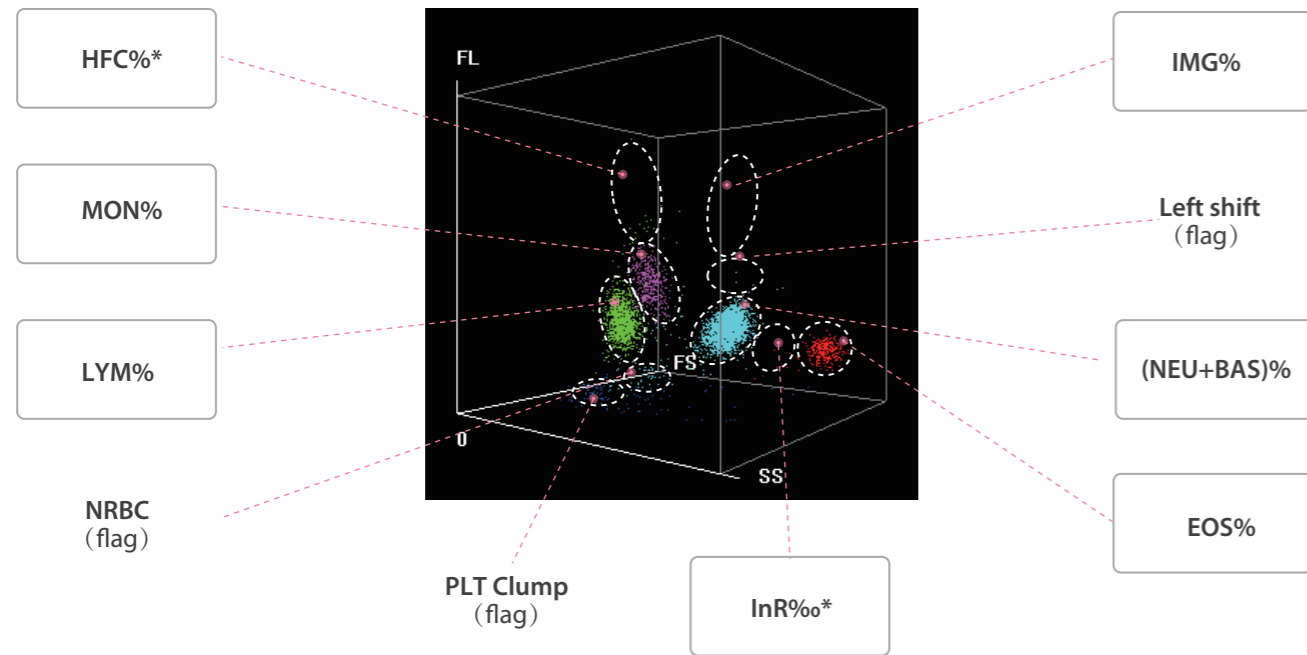

In DIFF scattergram, BC-6000 not only gives WBC 6-part differential results (with immature granulocyte), but also brings research parameters such as HFC (Blast & Atypical Lymphocyte), InR (information about malaria) and flags for Band, NRBC, PLT clump and Atypical Lymphocyte.


HFC\*(#, %) parameters represent high population of fluorescent cell, such as Blasts and Atypical Lymphocytes.

IMG(#, %) parameters provide information about immature granulocytes, including Promyelocytes, Myelocytes, Metamyelocytes, Immature Eosinophils and Immature Basophils.

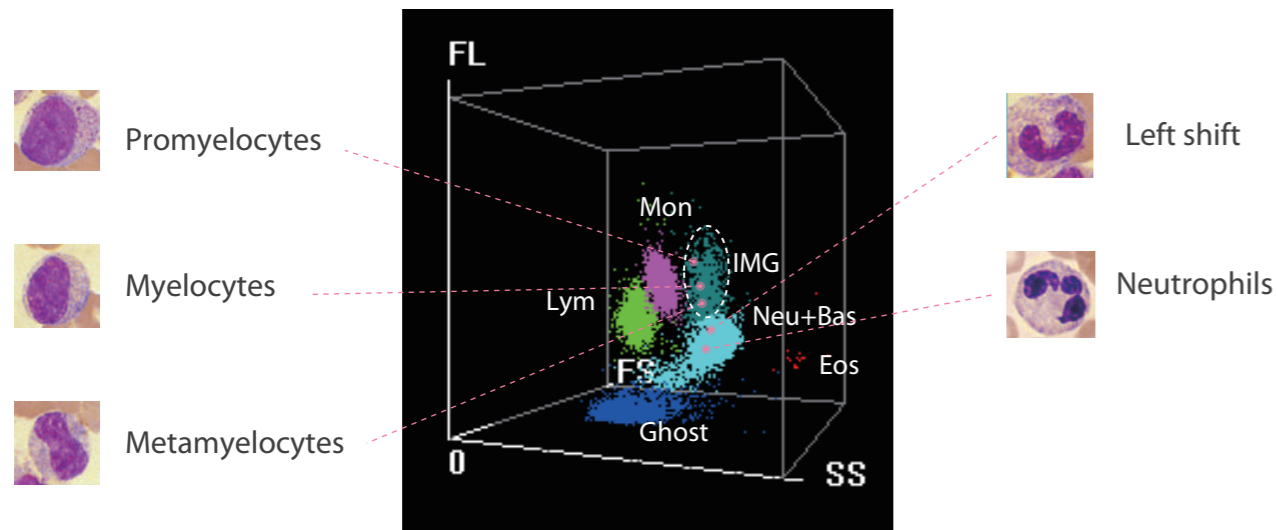

\*For research use only

# Malaria screening

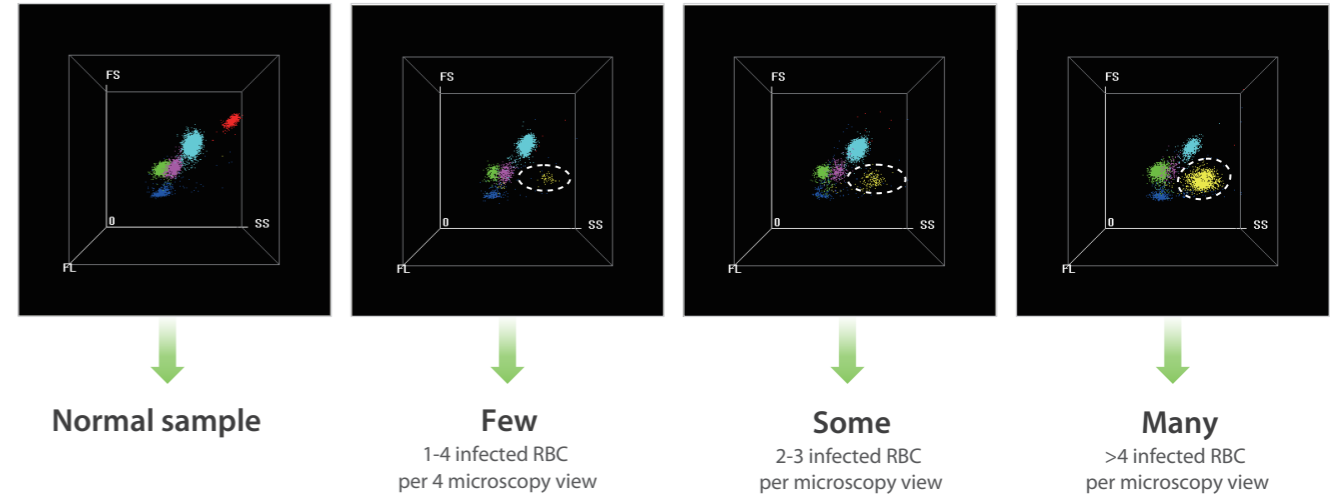


*Note: The yellow scatters are just for highlight.*

BC-6000 provides a dedicated flag called "infected RBC?", and "InR\*(#,%)" parameters to represent the number and ratio of the infected red blood cells in the sample respectively. BC-6000 users can obtain information about the possible presence of plasmodium parasite, the causative agent of malaria infection.

With the rising number of red blood cells with malaria parasites, the number of dots in the "InR" area increases proportionately. This creates the possibility to not only screen but also judge the severity of malaria infection.

# Interference prevention

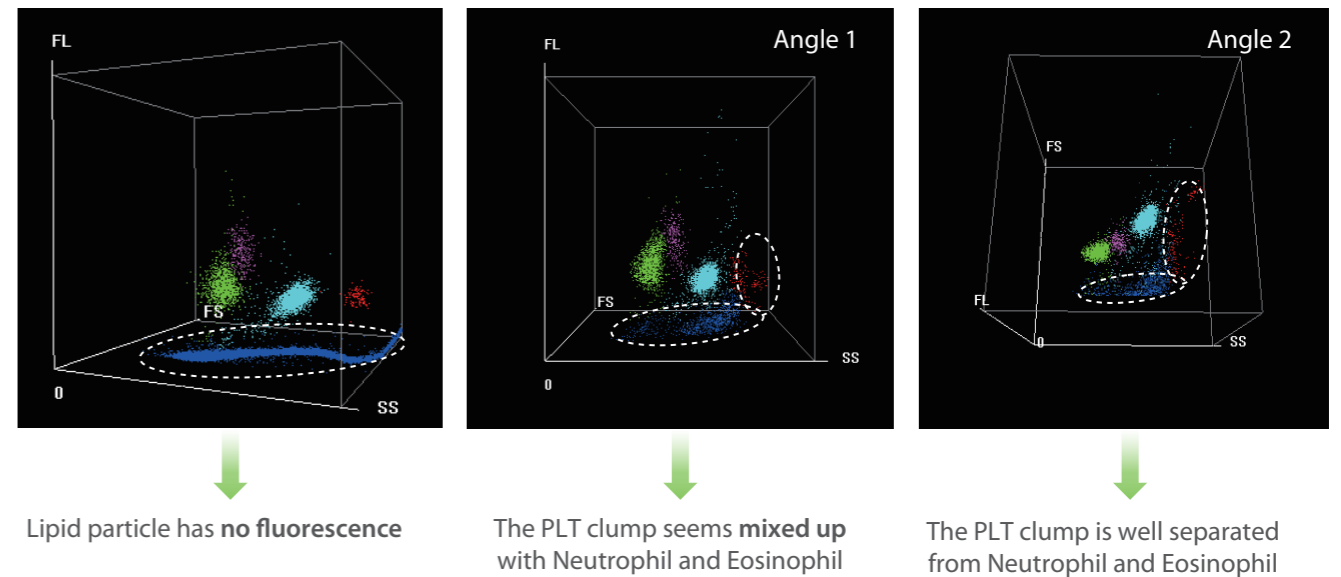


In DIFF scattergram, WBCs are dyed, but not lipid particles, by fluorescence, which prevents interference and ensures more accurate WBC results.

With information obtained through the 3D analysis, PLT clumps are well separated from each cluster of WBCs.



# WNB Channel

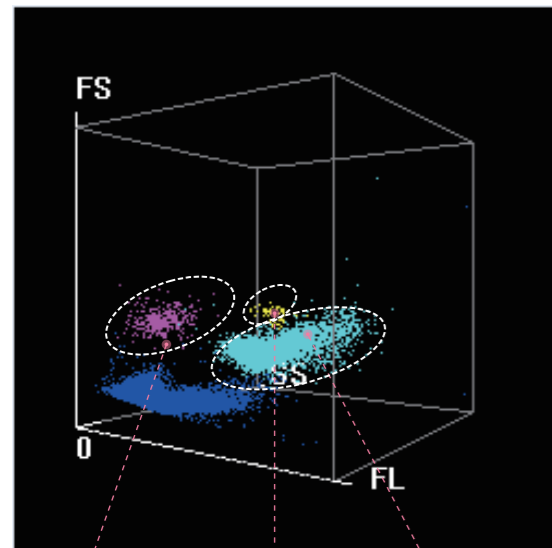

In WNB scattergram, BC-6000 provides NRBC, Basophils and WBC-N\* results. It means that the actual number of NRBCs can be measured in routine CBC, if they are present in the sample. Basophils are counted in this counting channel with NRBC results.

Basophil and NRBC results are generated on BC-6000 without extra reagent or cost.

NRBC Baso% WBC-N\*

*\*For research use only*

## NRBC results in every CBC




Automatic correction WBC counting, make sure neonatal counting correctly



Diagnosis for hemolytic anemia




Monitoring of hematopoietic diseases

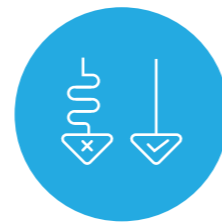

Reduce the ratio of review

NRBCs do not usually exist in the peripheral blood except that of newborn children. Detection of NRBCs is essential in diagnosing and monitoring the hematopoietic diseases.

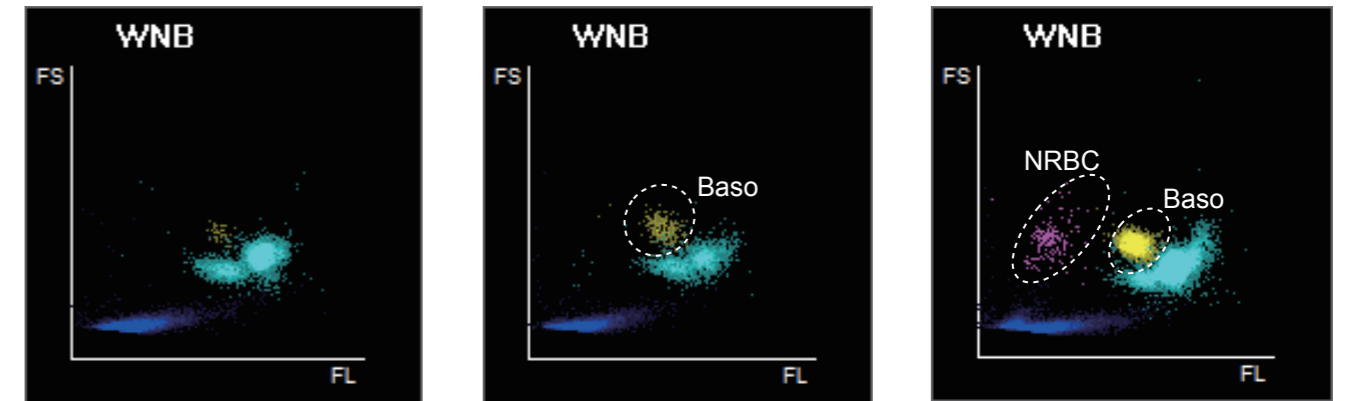


Normal sample

High Baso sample

High Baso & NRBC sample

BC-6000 provides accurate results on samples even with high level of Basophils and NRBCs.


## Body fluid

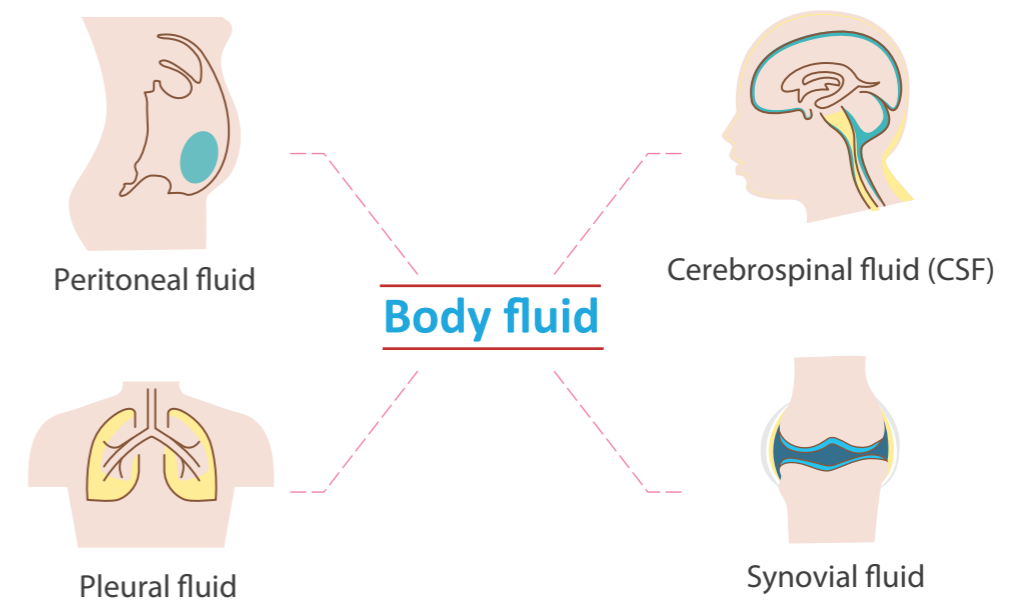

Besides blood specimen, BC-6000 also has body fluid test function without requiring dedicated reagent. The various types of body fluids include Peritoneal fluid, Pleural fluid, Cerebrospinal fluid (CSF) and Synovial fluid.

## Applicable to variety of tubes

To cater to customer's diversified needs, different types of blood collection tubes can be used on BC-6000, including regular whole blood vacuum tube, capillary blood microtainer tube and Sarstedt tube.

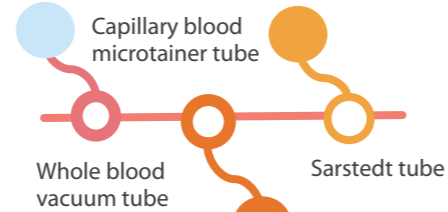

Operation

Cost



## Less testing time

BC-6000 can load up to 50 samples at a time and offers a throughput of up to 110 tests per hour.

## More Intuitive labXpert software

LabXpert is a standard configuration of BC-6000 for professional data analysis. The labXpert software optimizes functions to simplify your workflow for data analysis including improving re-exam efficiency, auto-validation for normal samples; it also provides more intuitive interface for you to review and validate pathological samples.




## Low sample volume

BC-6000 requires less sample volume as well as reagent consumption. For a CBC+DIFF test with NRBC result, BC-6000 only requires 80µL of whole blood and 35 µL of capillary blood.



## Easy maintenance

The only maintenance for end user is daily shut down by probe cleanser or probe cleanser cleaning once per day (if not shut down). The "auto-protect" program reminds operators when maintenance is needed (if not shut down).

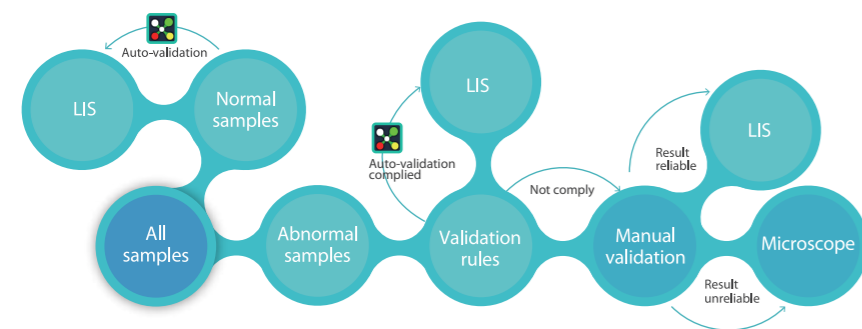

## Automatic Rerun & Reflex

Should the sample results trigger the criteria, the autoloader of BC-6000 can return the sample racks for an automatic rerun or reflex check.

# BS-240Pro

## Chemistry Analyzer

### Technical Specifications

#### System function

Automatic, Discrete, Random Access, Bench-top  
STAT sample priority

Throughput: Constant 240 photometric tests per hour, up to 400 T/H with ISE

Measuring principles: Absorbance photometry, turbidimetry, ion selective electrode technology

Methodology: End-point, Fixed-time, Kinetic, optional ISE, Single/Double reagent chemistries, Mono-chromatic / bi-chromatic

Original system pack reagent ready to use

Close system and open system is optional

#### Reagent/Sample Handling

Reagent/Sample tray: 50 to 100 positions for reagents and 50 to 100 positions for samples in 24-hour refrigerated compartment (2~12°C)

Reagent volume: R1: 100~200µL, step by 0.5µL

R2: 10~200µL, step by 0.5µL

Sample volume: 2~35µL, step by 0.1µL

Reagent/Sample probe: Liquid level detection, horizontal and vertical collision protection, inventory checking, reagent pre-warming, optional clog detection

Probe cleaning: Automatic washing for interior and exterior  
Carry over < 0.05%

Automatic sample dilution: Pre-dilution and post-dilution

Mixing Unit: Independent mixing bar

#### Built-in Bar Code Reader (Optional)

Used for sample and reagent programming

Be applicable to various bar code systems of Codabar, ITF (Interleaved Two of Five), code128, code39, UPC/EAN, Code93

Capable to communicate with LIS in bi-directional mode

#### Reaction System

Reaction tray: 80 reusable cuvettes

Reaction volume: 100~360µL

Reaction temperature: 37°C ± 0.1°C

Cuvette Washing: Washing station with pre-warmed detergent and de-ionized water

#### ISE Module (optional)

Measuring K<sup>+</sup>, Na<sup>+</sup>, Cl<sup>-</sup>

#### Optical System

Light Source: Halogen-tungsten lamp

Wavelength: 12 wavelengths, 340nm, 380nm, 412nm, 450nm, 505nm, 546nm, 570nm, 605nm, 660nm, 700nm, 740nm, 800nm

Absorption range: 0~3.5Abs, resolution 0.0001Abs

Stray Light: 4.9Abs

#### Control and Calibration

Calibration modes: K factor, Linear (two points and multi-points), Logit-Log 4P, Logit-Log 5P, Spline, Exponential, Polynomial, Parabola, Logit-Log3P, Broken line

Control Rules: One key calibrator import function  
Westgard multi-rule, Levey-Jennings, Cumulative sum check, Twin plot

#### Operation Unit

Operation system: Windows 10

Interface: RS-232

#### Working Conditions

Power Supply: 200~240V, 50/60Hz, ≤1300VA or 100~130V, 60Hz, ≤1300VA

Dimension: 860 mm (length) × 660 mm (depth) × 550 mm (height)

Weight: 115 kg

Water Consumption: ≤6.5 L/H



## BS-240Pro

### Chemistry Analyzer

# Compact yet Robust

Mindray Building, Keji 12th Road South,  
High-tech Industrial Park, Nanshan, Shenzhen 518057, P.R. China  
Tel: +86 755 8188 8998 Fax: +86 755 26582680  
E-mail: intl-market@mindray.com www.mindray.com

**mindray** and **BS-240Pro** are registered trademarks or trademarks owned by Shenzhen Mindray Bio-medical Electronics Co., Ltd.  
© 2015 Shenzhen Mindray Bio-Medical Electronics Co., Ltd. All rights reserved. Specifications subject to changes without prior notice.  
P/N: ENG-BS-240Pro-210285X8P-20170703

**mindray**

**mindray**  
healthcare within reach



## Large and flexible capacity

Up to 100 sample positions  
Up to 100 reagent positions  
(50 fixed + 50 interchangeable)



## Gratings photometer

100 $\mu$ L minimum reaction volume



## A Whole New Generation

with constant throughput  
of 240 photometric T/H

## Considerate design

More space to load samples



## HbA1c smart-sampling function

Automatic preparation of hemolysate

# BS-240Pro


Chemistry Analyzer



Waterfall probe cleaning



Intelligent probe with optional clog detection

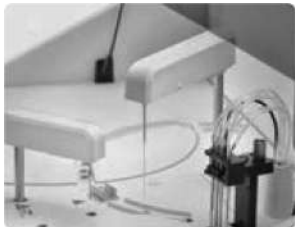

Constant throughput



Independent mixing bar



Optimized washing station



Built-in barcode reader



Optional ISE module easy to access



Intuitive software with more functionalities




## Complete traceability process

Complete calibration hierarchy and traceability chain are based on ISO standard (EN/ISO17511) from reference system to routine measurement system.

### Traceability chain of Mindray measurement system (Glu)

Traceability Material Calibration / Value Assignment Procedure Implementation Uncertainty U<sub>cj</sub>

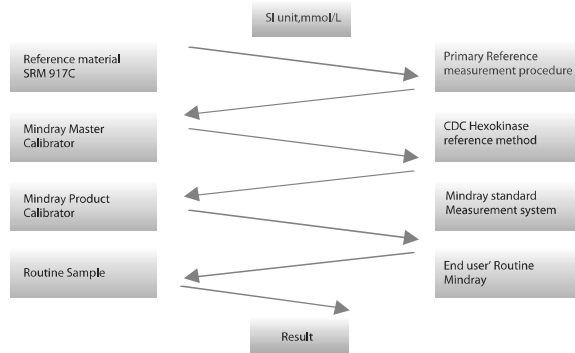


## External quality assurance for reference measurement

Mindray participates in RELA (External quality control for reference laboratory).

## EQA for Mindray Reference laboratory——RELA

Mindray reference laboratory has passed RELA for 6 consecutive years.

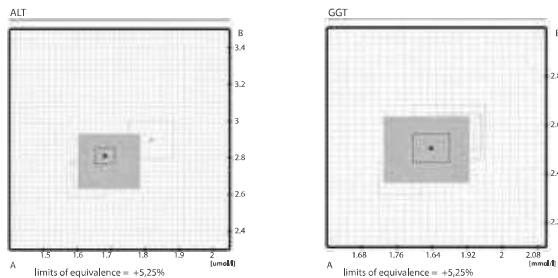

More RELA results please refer to: [www.dgkl-rfb.de/81](http://www.dgkl-rfb.de/81)



## Reagent menu

### Hepatic Panel

Alanine Aminotransferase (ALT)  
Aspartate Aminotransferase (AST)  
Alkaline Phosphatase (ALP)  
 $\gamma$ -GlutamylTransferase ( $\gamma$ -GT)  
Direct Bilirubin (D-Bil) DSA Method  
Direct Bilirubin (D-Bil)VOX Method  
Total Bilirubin (T-Bil) DSA Method  
Total Bilirubin (T-Bil)VOX Method  
Total Protein (TP)  
Albumin (ALB)  
Total Bile Acids (TBA)  
Prealbumin (PA)  
Cholinesterase (CHE)  
 $\alpha$ -L-fucosidase (AFU)  
5'-nucleotidase (5'-NT)

### Renal Panel

Urea (UREA)  
Creatinine (CREA)Modified Jaffé Method  
Creatinine (CREA)Sarcosine Oxidase Method  
Uric Acid (UA)  
Carbon dioxide (CO2)  
Microalbumin(MALB)  
 $\beta$ 2-Microglobulin ( $\beta$ 2-MG)  
Cystatin C (CysC)  
Retinol binding protein( RBP)

### Immune Panel

Immunoglobulin A (IgA)  
Immunoglobulin G (IgG)  
Immunoglobulin M (IgM)  
Complement C3 (C3)  
Complement C4 (C4)

### Diabetes Panel

Glucose (Glu) GOD-POD Method  
Glucose (Glu) HK Method  
Hemoglobin A1c (HbA1c)  
Fructosamine (FUN)  
 $\beta$ -Hydroxybutyrate( $\beta$ -HB)

### Cardiac panel

Creatine Kinase (CK)  
Creatine Kinase-MB (CK-MB)  
Lactate Dehydrogenase (LDH)  
 $\alpha$ -Hydroxybutyrate Dehydrogenase( $\alpha$ -HBDH)  
High sensitive C-reaction protein( HS-CRP)

### Inorganic & Anemia

Iron (Fe)  
Ferritin (FER)  
Transferrin (TRF)  
Calcium (Ca)  
Magnesium (Mg)  
Phosphate Inorganic (P)  
Unsaturated iron binding capacity (UIBC)  
Glucose-6-phosphate dehydrogenase (G6PD)

### Lipid Panel

Total Cholesterol (TC)  
Triglycerides (TG)  
HDL-Cholesterol (HDL-C)  
LDL-Cholesterol (LDL-C)  
Apolipoprotein A1 (ApoA1)  
Apolipoprotein B (ApoB)  
Lipoprotein(a) [Lp(a)]

### Rheumatism Panel

C-reactive protein (CRP)  
Rheumatoid Factor (RF)  
Antibodies Against Streptolysin O (ASO)

### Lung Panel

Adenosine Deaminase (ADA)  
Angiotensin Converting Enzyme(ACE)

### Pancreatitis Panel

$\alpha$ -Amylase ( $\alpha$ -AMY)  
Lipase (LIP)

## TECHNICAL SPECIFICATIONS

<b>Principle</b>	Ion selective Electrode ( ISE )
<b>Sample Types</b>	Serum, Plasma, Whole Blood, Urine, CSF
<b>Parameters</b>	Na <sup>+</sup> , K <sup>+</sup> , Cl <sup>-</sup> , Ca <sup>++</sup> & Li <sup>+</sup>
<b>Sample volume</b>	≤ 90μL
<b>Analysis Time</b>	≤ 60 secs
<b>Throughput</b>	≤ 60 samples per hour
<b>Calibration</b>	Automatic Calibration : <ul style="list-style-type: none"> <li>• Two Point Calibration after 4 hrs or just before Sample analysis</li> <li>• One Point Calibration with each sample analysis</li> </ul>
<b>Normal Ranges</b>	User defined normal ranges for all sample types with flagging for abnormal samples
<b>Data Management</b>	Calibration History, Sample History and control history of last 30 days
<b>Quality Control</b>	Tri-Level Quality Control, with L J Plots of last 30 days
<b>Display</b>	5" color LCD with Touchscreen Display, Resolution: 480 x 272 pixels
<b>Printer</b>	Built-in 2" ( 48 mm ) thermal Printer, 42 characters wide
<b>External Interface</b>	USB ( Keyboard, Mouse )
<b>Operating Temp.</b>	15 to 35°C
<b>Relative Humidity</b>	5 to 85% non-condensing
<b>Power Supply</b>	100 to 240VAC 50/60 Hz
<b>Dimensions</b>	395mm x 260mm x 375mm
<b>Weight   Sample Storage</b>	7.5 Kg   100,000

\* Optional Barcode reader available.

### ISE - Tritrol

- Ready to use assayed control
- 3 different levels of concentrations - Low, Normal & High
- 1 year closed vial stability

#### Pack size Available : 13mL

- Aqueous control - Low, Normal & High ( Each of 3mL )
- Urine control - Normal & High ( Each of 1mL )
- IFS ( 1 X 2mL )



CE

ERMA INC.  
Yushima 2-31-6,  
Bunkyo-ku,  
Tokyo, Japan

Tel. : 81-3-3818-6281  
Fax. : 81-3-3818-6934  
Email : international@erma.co.jp  
Web : www.erma.jp

\*Specifications are subject to change without notice.



ACCURATE | RELIABLE | FRIENDLY INTERFACE



# EL-120

## ELECTROLYTE ANALYZER

- Na<sup>+</sup>, K<sup>+</sup>, Cl<sup>-</sup>, Ca<sup>++</sup>, Li<sup>+</sup>
- Low Sample Volume
- Fast Sample Processing
- Low Cost per Test

www.erma.jp

Erma EL - 120 is ideal for medium to high volume Laboratories and Hospitals. It is user friendly & reliable. It also allows measurement of specific ions; accurately and precisely. This supports fast diagnosis.

## MAIN MENU

- Large Graphical Display affords convenience of operation and minimizes errors while entering data
- User friendly software with interactive Graphical user interface
- Shortcut Menu Keys for Easy Operation



## ANALYSIS

- Aspiration system: Internal peristaltic Stepper Motor Pump for accurate and precise aspiration of analytes
- Automatic flushing after aspiration of every sample to eliminate carryover
- Precise & accurate result



## ELECTRODES

- Electrodes are based on Ion Selective principle
- Electrodes can be refilled.
- Long life and assured stability



## CALIBRATION

- Two Point Calibration after 4 hrs/ 8 hrs or 12 hrs. Programmable as per user
- One Point Calibration with each Sample Analysis
- Displays status of calibration with date and time for each parameter
- Also displays the date and time of next calibration to be performed / due

## RUN QUALITY CONTROL

- Levy Jennings plot of daily and monthly QC data to assess the performance
- Helpful for laboratories incase of accreditation and quality assessment

## ISEPak

- Convenient easy to place ISEpak which contains reagents required for all parameters
- Integrated reagent Pack for testing of all parameters
- In-Built reagent inventory management software

### Measuring Range

PARAMETER	SERUM / PLASMA BLOOD / AQUEOUS
Na <sup>+</sup>	40 - 205
K <sup>+</sup>	1.5 - 15
Cl <sup>-</sup>	50 - 200
Ca <sup>++</sup>	0.2 - 5.0
Li <sup>+</sup>	0.1 - 6.0

Control Result  
 Date : 2018-11-30 15:04:08  
 Lot No : 5  
 Control : LEVEL 2  
 Sample Type : Aqueous

Param	mmol/L	Target
Na <sup>+</sup>	141.3	140 (4.0)
K <sup>+</sup>	4.33	4.25 (0.3)
Cl <sup>-</sup>	104.7	105 (4.0)
Ca <sup>++</sup>	1.31	1.25 (0.2)





# MAGLUMI 800

Chemiluminescence Immunoassay (CLIA) System



# OUTSTANDING TECHNOLOGY POWER OF MAGLUMI

## Chemiluminescence Immunoassay (CLIA) System

CLIA uses two important technologies, one is labeling technology which determines reaction mode; and the other is separation technology which determines the sensitivity, accuracy and precision of the reagents.

### Labeling Technology

Two types of labeling technologies are commonly used. One is enzyme label and the other is non-enzyme small molecule label. Enzyme label reagents are not so stable and are easily effected by the change of storage conditions. Maglumi system applies ABEI labels. ABEI is a non-enzyme small molecule with special molecular formula to enhance stability in acid and alkaline solution. As fast chemiluminescence, ABEI chemical reaction with sodium hydroxide (NaOH) and Hyperoxide ( $H_2O_2$ ) finishes the process in 3 seconds.

### Separation Technology

Maglumi uses Magnetic Microbeads. As separation technology, it has been widely used in the field of CLIA. Compared with traditional separation technology, it has the following advantages:

- Shortening the reaction time by enlarging the reaction area of antigens and antibodies.
- Enhancing the sensitivity by better and faster capture of antigens and antibodies.
- Reducing inter or intra-assay discrepancies significantly by mixing the reagents thoroughly in a liquid separation platform.
- Enhancing the accuracy by absorbing antigens and antibodies through chemical reaction.



# ADVANTAGE OF MAGLUMI 800




- For small and middle size labs or hospitals
- Simple, Smart, Humanization designed model
- More than enough test menu to meet the different needs and it is enlarging according to your requirement
- Free Quality Control and calibrator helps to reduce the Cost Per Test
- Quality guarantee by the Third party quality control and External Quality Assessment (EQA)
- Key technology ABEI and Nano magnetic microbead enhances stability and sensitivity of MAGLUMI reagents




# MAGLUMI 800

- Integrated design and humanized appearance
- Simple and smart operation
- Compact size and weight

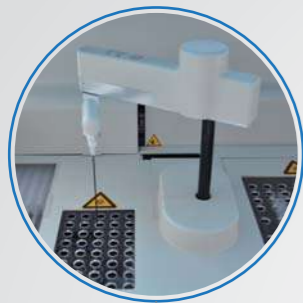


## Pipettor

- Aspirate reagent/sample with high speed
- Titanium needle for crushproof
- Clot detection & liquid level detection
- Coated with TEFLON to prevent carry-over
- Inner & external washing

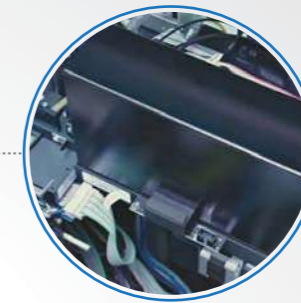


## Chamber

- High-sensitivity and low-noise PMT
- Cuvettes detection and anti-overflow

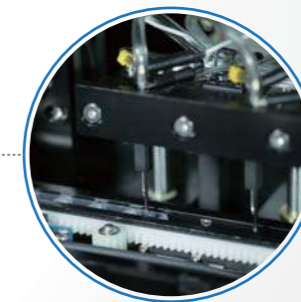

## Washer

- Strong magnetic field helps better separation
- Three pairs of washing needles with anti-overflow function

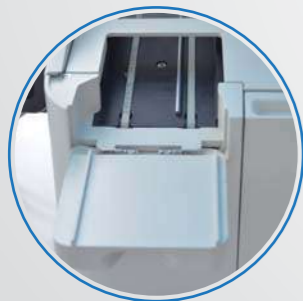


## Automated cuvette loader

- Up to 40 modules of 240 tests
- Continuous loading




## Sample area

- Up to 40 sample on board with barcode
- Continuous loading & cooling function

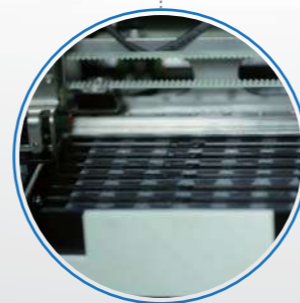

## Incubator

- Incubate 13 slots (78 tests at the same time)
- The temperature is  $36.8\text{ C} \pm 0.5\text{ C}$

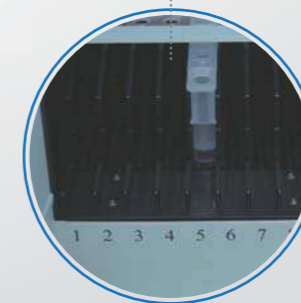

## Reagent area

- Up to 9 reagents on board with RFID
- RFID label storing reagent information
- Easy and fast for adding and updating
- RFID tag, high reliability

# OPERATION SOFTWARE



## Comprehensive software

- User-friendly Interface
- Real-time status monitoring for each test
- Monitoring reagent and consumable status
- Intelligent alarm function



## Maintenance Guide

- Daily, weekly and monthly maintenance guide for analytical unit
- Ensure performance reliability and reduce unnecessary service calls



## Quality control

- Westgard rules and Levey-Jennings chart for both internal and external quality control



## Test Summary Function

- Test summary including system test, calibration, QC, statistics of samples, valid tests and rerun tests
- Search and review test information conveniently



## Pre or post defined dilution Function

- Pre or post defined dilution flexibly if high concentration is met



## Connect to LIS ( bidirection )

- Bidirection, Lis can connect by online and ASTM

# TEST MENU

## Tumor Markers

Ferritin  
AFP  
CEA  
Total PSA  
f-PSA  
CA 125  
CA 15-3  
CA 19-9  
PAP  
CA 50  
CYFRA 21-1  
CA 242  
CA 72-4  
NSE  
S-100  
SCCA  
TPA-snibe  
Pepsinogen I  
Pepsinogen II  
Gastrin-17  
H.pylori IgG

## Bone Metabolism

Intact PTH  
Calcitonin  
Osteocalcin  
25-OH Vitamin D

## Glyco Metabolism

C-Peptide  
Insulin  
ICA  
IAA  
Proinsulin  
GAD 65  
IA-2

## Fertility

FSH  
LH  
HCG/ $\beta$ -HCG  
PRL  
Estradiol  
free Estriol  
Progesterone  
Testosterone  
free Testosterone  
DHEA-S  
**\*AMH**  
17-OH Progesterone

## Thyroid

TSH  
T4  
T3  
FT4  
FT3  
TG  
TGA  
TRAb  
TMA  
Anti-TPO  
Rev T3

## Prenatal Screening

cAFP  
free  $\beta$ -HCG  
PAPP-A

## Inflammation Monitoring

CRP  
PCT

## Kidney Function

$\beta$ 2-MG  
Albumin

## Hepatic Fibrosis

HA  
PIIIP N-P  
C IV  
Laminin  
Cholyglycine

## TORCH

Toxo IgG  
Toxo IgM  
Rubella IgG  
Rubella IgM  
CMV IgG  
CMV IgM  
HSV-1/2 IgG  
HSV-2 IgG  
HSV-1/2 IgM

## EBV

EBV EA IgG  
EBV EA IgA  
EBV VCA IgG  
EBV VCA IgM  
EBV VCA IgA  
EBV NA IgG

## Immunoglobulin

IgM  
IgA  
IgE  
IgG

## Anemia

Vitamin B12  
Ferritin  
FA

## Drug Monitoring

Cyclosporine A  
Tacrolimus, FK 506  
Digoxin

## Cardiac

CK-MB  
Troponin I  
Myoglobin  
NT-proBNP  
Aldosterone  
Angiotensin I  
Angiotensin II  
D-Dimer  
**\*Direct Renin**  
LP-PLA2  
hs-cTnI

## Infectious Disease

HBsAg  
Anti-HBs  
HBeAg  
Anti-HBe  
Anti-HBc  
Anti-HCv  
Syphilis  
Chagas  
HTLV I/II  
Anti-HAV  
HAV IgM  
HIV A b  
HIV p24 Ag  
HIV Ab/Ag combi

## Others

GH  
IGF-I  
Cortisol  
ACTH

**\*In Red Available Soon**

# MAGLUMI 800

## Chemiluminescence Immunoassay (CLIA) System

### Technical Specifications

#### Reagent Feature

- Flash chemiluminescence label-ABEI,with high sensitivity and long stability
- Nano magnetic microbeads separation,fast and efficient

#### Main Features

- Throughput: up to 180 tests/hour
- 24 hours ready to use
- Time to first result: 17 minutes

#### Modes of Operation

- Random,Batch and STAT

#### Sample Loading

- Up to 40 sample tubes
- Continuous loading,STAT available
- Barcode reader recognition or analyzer auto numbered
- LIS connection,auto read sample info
- Refrigerated sample area with independent power supply

#### Reagent Loading

- 9 reagents on board
- Continuous loading
- RFID reading all info of reagent
- Refrigerated reagent area

#### Reagent Integral

- Ready-to-use, no pretreatment required
- Including calibrators
- RFID tag storing all info of reagent
- RFID tag with built-in master curve
- 2-point calibration to adjust master curve
- Calibration stability: Max 4 weeks

#### Other Features

- Clot detection
- Liquid level detection
- Auto and optional ratio for high sample dilution
- Constant 37°C incubation

#### Operation System

- Windows 7
- DELL compatible (Dual-core CPU)
- Color touch screen

#### Interconnection

- Bi-communication with LIS via ASTM protocol

#### Dimensions Weight

- Input:AC 110-230V,50/60 Hz,Max 630W
- Dimension and Weight:102\*72\*56cm 73kg



Shenzhen New Industries Biomedical Engineering Co., Ltd.(Snibe)

21st Floor, Block A, Building 1, Shenzhen Software Industry Base,  
No. 1008, Keyuan Road, Nanshan District, Shenzhen 518000 CHINA

Tel: +86 755 26501514

Fax: +86 755 26654850

Email: sales@snibe.com

Web: www.snibe.com



DISTRIBUTOR:

# MAGLUMI<sup>®</sup> X3


Fully-auto Chemiluminescence Immunoassay (CLIA) System



## Save Your Space without Compromise

Powerful Performance and Compact Design  
Compatible with Small and Medium-sized Hospitals and Labs



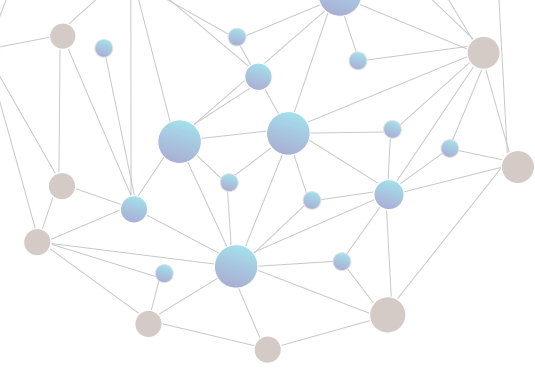


# MAGLUMI® X3

## Features & Benefits

- Small but powerful, the throughput is up to 200 tests/hour, and throughput per unit area is 294 T/h/m<sup>2</sup>.
- Compatible with all MAGLUMI® reagents, one of the broadest automated CLIA test menus in the world (166 parameters).
- The latest intelligent washing technology and bidirectional temperature control measurement guarantee accurate and reliable results. The comprehensive advanced design of MAGLUMI® X3 ensures excellent performance.
- Single reaction cup can avoid light pollution and increase cuvette utilization, its integrated packaging can avoid the stuck of the cuvette and destroying the cuvette wall.
- No-pause loading/unloading of reagents/samples without waiting or interrupting tests.

## Functional Modules

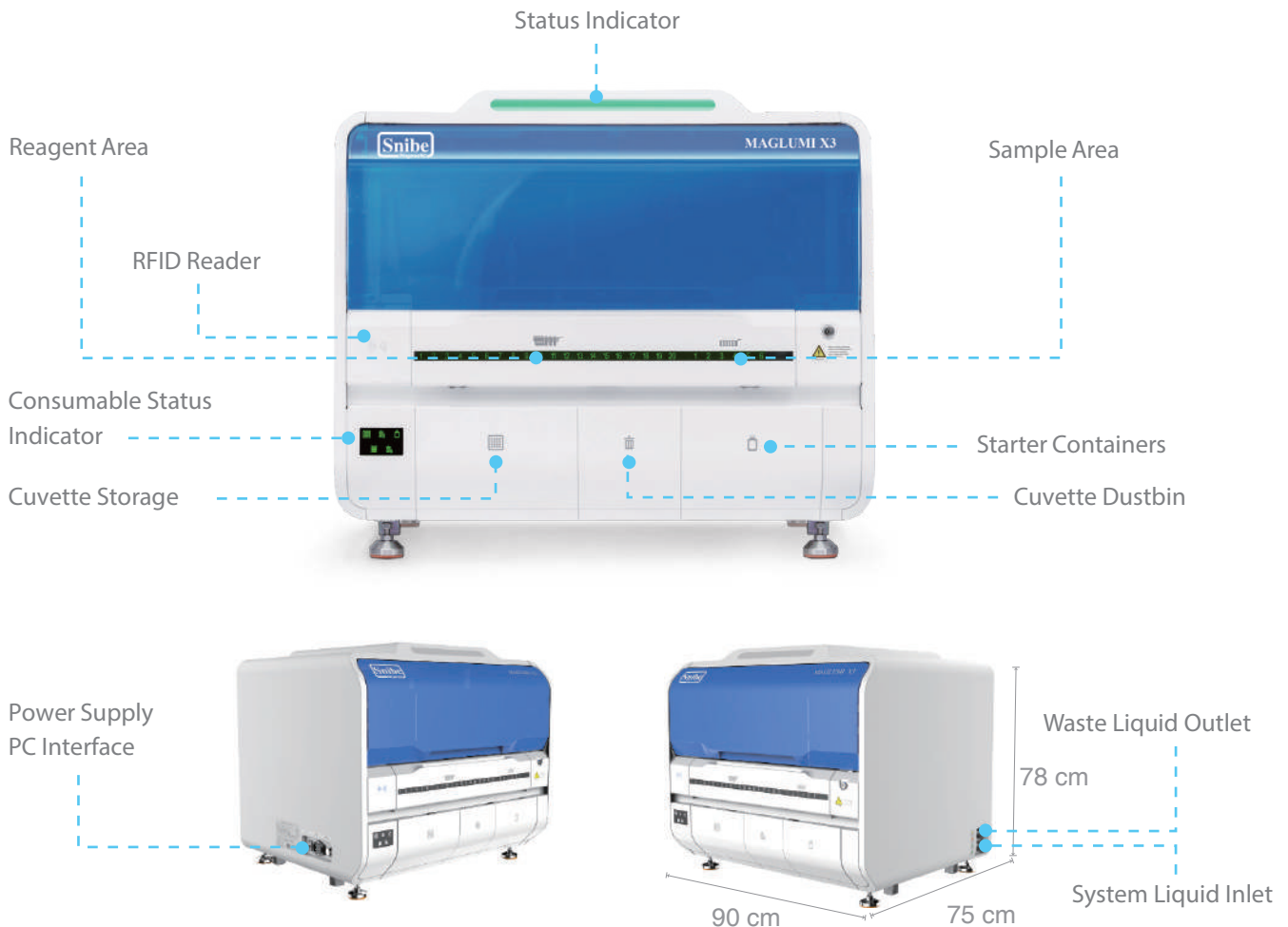




## Reagent and Sample Loading

- 72 sample positions with barcode reader
- Sample in each position can be defined as STAT
- 20 reagent positions with RFID reader, 24h refrigerated at 10°C ( $\pm 0.3^{\circ}\text{C}$ )
- No-pause loading/unloading of reagents/samples without waiting or interrupting tests
- Reagent and sample indicator light available, real-time status monitoring, no need to check on the computer monitor



## Single Reaction Cup with Integrated Packaging

- Single reaction cup can avoid light pollution and increase the cuvette utilization
- Single reaction cup with integrated packaging can avoid the stuck of the cuvette and destroying the cuvette wall

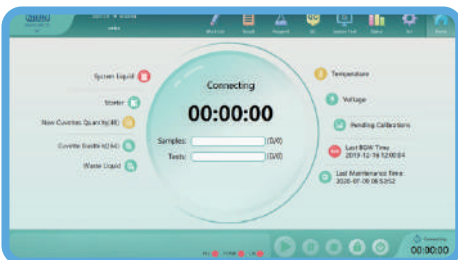

## Intelligent UI Interface

- New user-friendly UI interface (Support multi-language)
- Real-time display for monitoring each test procedure, reagent, and consumable quantities, it's convenient for the management
- Intelligent error recovery



## Convenient Management of Analyzer Status

- Intuitive indicator light of reagent, sample, and consumables, no need to focus on the monitor, the status can be known by looking at the analyzer from a distance



# All Balanced and Strong

## Accurate Pipetting Technology

- Single needle coated with TEFLON to avoid carry-over
- Independent washing unit for the needle, with different washing modes depending on different assay
- Pipetting detection technologies include liquid level, clot, and crash detection can ensure accurate and normal pipetting




## Accurate Incubation Technology

- Single reaction cup with five-sided heating at 37.0°C ( $\pm 0.3^\circ\text{C}$ ), ensuring accurate, quick and uniform incubation
- Incubating 80 cuvettes simultaneously
- Non-contact mixing with different modes before incubation, ensuring sufficient reaction



## High-Efficiency Washing Technology

- Magnetic separation washing with 4-step independent units and non-contact vortex mixing can reduce non-specific binding
- Different mixing modes depending on different assays, ensuring a sufficient washing to improve sensitivity

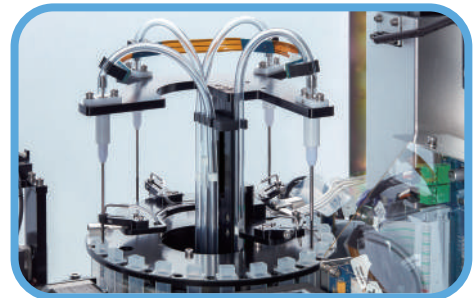

## Stable and Accurate Measuring Chamber


- ABEI-label exclusive proprietary technology
- PMT with heating and cooling function, ensuring the detection precision
- Independent measuring room, avoiding light cross-contamination






## Exclusive Intelligent Washing Technology

### Independent washing units for waste liquid needles

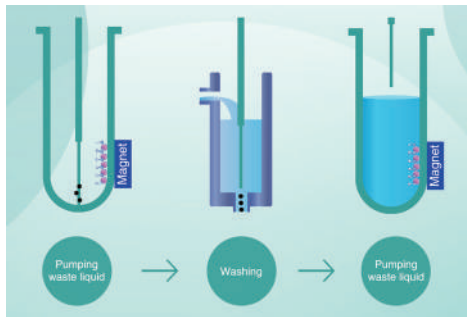




In the wash station, four waste liquid needles are equipped with four independent washing units to avoid the carry-over.

### Magnetic separation washing and non-contact vortex mixing

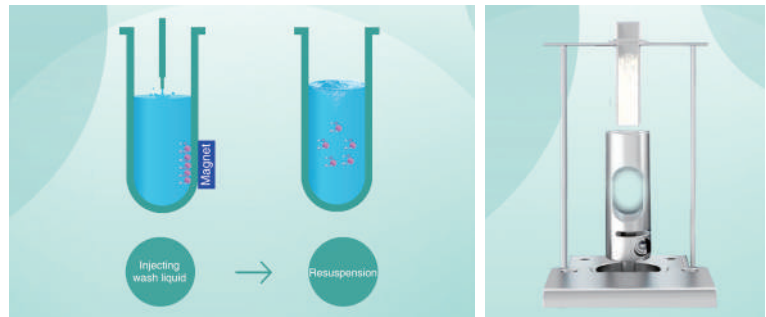

Magnetic separation washing and non-contact vortex mixing units, with different mixing modes depending on different assays, ensuring a sufficient washing to reduce the non-specific binding.

## Stable and Accurate Measuring Chamber

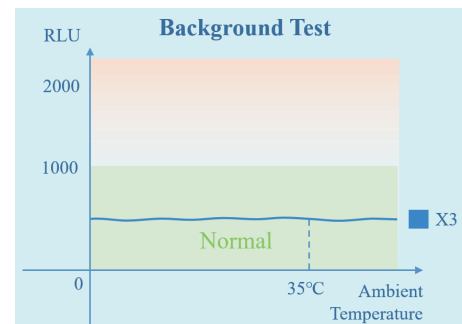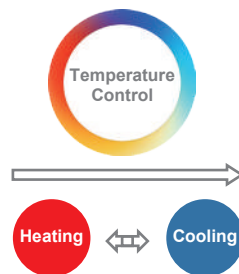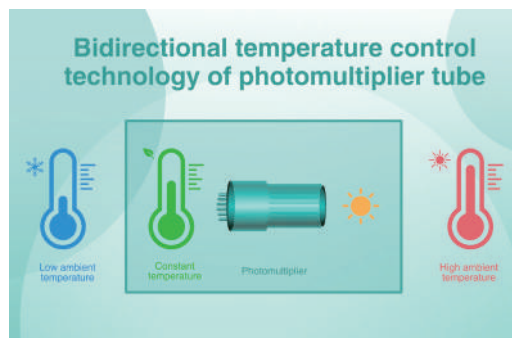

Automatic temperature control in the measuring chamber can keep a stable background value, not affected by the ambient temperature, ensuring the accuracy of measurement results.



## Top 10 Highlights of MAGLUMI® X3

- 1 One of the highest space efficiency CLIA analyzers combines one of the broadest automated CLIA test menus.
- 2 The throughput is up to 200 tests/hour (Throughput per unit area is 294 T/h/m<sup>2</sup>). No worry about the lack of laboratory space.
- 3 **MAGLUMI®** reagent kits (166 parameters) and consumables are compatible with **MAGLUMI® X3**.
- 4 The comprehensive advanced design of **MAGLUMI® X3** ensures excellent performance.
- 5 Single reaction cup can avoid light pollution and increase cuvette utilization, its integrated packaging can avoid the stuck of the cuvette and destroying the cuvette wall.
- 6 No-pause loading/unloading of reagents/samples without waiting or interrupting tests.
- 7 Pipetting detection technologies include liquid level, clot, and crash detection can ensure accurate and normal pipetting.
- 8 The latest intelligent washing technology and bidirectional temperature control measurement guarantee accurate and reliable results.
- 9 Intelligent UI interface, user-friendly, easy, and convenient operation.
- 10 Intuitive indicator light make users no need to check reagents and consumables frequently.

## MAGLUMI® X3 Reagent and Consumables



Reagent Kit



Reaction Module



Starter 1+2



Light Check



Wash Concentrate



System Tubing Cleaning Solution

# Broad CLIA Test Menu with 166 Assays



## Thyroid

TSH (3rd Generation)  
T4  
T3  
FT4  
FT3  
Tg (Thyroglobulin)  
TGA (Anti-Tg)  
Anti-TPO  
TRAb  
TMA  
Rev T3  
\*T-Uptake

## Hepatic Fibrosis

HA  
PIIIP N-P  
C IV  
Laminin  
Cholyglycine

## TORCH

Toxo IgG  
Toxo IgM  
Rubella IgG  
Rubella IgM  
CMV IgG  
CMV IgM  
HSV-1/2 IgG  
HSV-1/2 IgM  
HSV-2 IgG  
\*HSV-2 IgM  
\*HSV-1 IgG  
\*HSV-1 IgM

## EBV

EBV EA IgG  
EBV EA IgA  
EBV VCA IgG  
EBV VCA IgM  
EBV VCA IgA  
EBV NA IgG  
EBV NA IgA

## Inflammation Monitoring

CRP (High Sensitive)  
PCT (Procalcitonin)  
IL-6 (Interleukin 6)  
SAA (Serum Amyloid A)

## Fertility

FSH  
LH  
HCG/ $\beta$ -HCG  
PRL (Prolactin)  
Estradiol  
Testosterone  
free Testosterone  
DHEA-S  
Progesterone  
free Estriol  
17-OH Progesterone  
AMH  
SHBG  
Androstenedione  
\*PIGF  
\*sFlt-1

## Hypertension

Direct Renin  
Aldosterone  
Angiotensin I  
Angiotensin II  
Cortisol  
ACTH

## Autoimmune

Anti-CCP  
Anti-dsDNA IgG  
ANA Screen  
ENA Screen  
Anti-Sm IgG  
Anti-Rib-P IgG  
Anti-Scl-70 IgG  
Anti-Centromeres IgG  
Anti-Jo-1 IgG  
Anti-M2-3E IgG  
Anti-Histones IgG  
Anti-nRNP/Sm IgG  
Anti-SS-B IgG  
Anti-SS-A IgG  
TGA(Anti-Tg)  
Anti-TPO  
TRAb  
TMA  
ICA  
IAA(Anti Insulin)  
GAD 65  
Anti-IA2  
\*Anti-Cardiolipin IgG  
\*Anti-Cardiolipin IgM  
\*Anti-MPO

## Immunoglobulin

IgM  
IgA  
IgE  
IgG

## Tumor Markers

AFP  
CEA  
Total PSA  
f-PSA  
CA 125  
CA 15-3  
CA 19-9  
PAP  
CA 50  
CYFRA 21-1  
CA 242  
CA 72-4  
NSE  
S-100  
SCCA  
TPA-snibe  
ProGRP  
HE4  
HER-2  
PIVKA-II

## Infectious Disease

HBsAg  
Anti-HBs  
HBeAg  
Anti-HBe  
Anti-HBc  
Anti-HCV  
Syphilis  
Anti-HAV  
HAV IgM  
HIV Ab/Ag Combi  
Chagas  
HTLV I+II  
H.pylori IgG  
H.pylori IgA  
H.pylori IgM  
2019-nCoV IgG  
2019-nCoV IgM  
SARS-CoV-2 S-RBD IgG  
SARS-CoV-2 Neutralizing Antibody  
SARS-CoV-2 Ag  
\*Anti-HBc IgM

## Drug Monitoring

Digoxin  
CSA (Cyclosporine A)  
FK 506 (Tacrolimus)

## Kidney Function

$\beta_2$ -MG  
Albumin  
\*NGAL

## Cardiac

CK-MB  
Troponin I  
Myoglobin  
hs-cTnI  
H-FABP  
NT-proBNP  
BNP  
D-Dimer  
Lp-PLA2  
\*MPO

## Metabolism

Pepsinogen I  
Pepsinogen II  
Gastrin-17  
GH (hGH)  
IGF-I  
IGFBP-3

## Prenatal Screening

AFP(Prenatal Screening)  
Free  $\beta$ -HCG  
PAPP-A  
free Estriol

## Anemia

Vitamin B12  
Ferritin  
Folate (FA)  
\*RBC Folate  
\*EPO

## Glyco Metabolism

C-Peptide  
Insulin  
GAD 65  
Anti-IA2  
ICA  
IAA (Anti Insulin)  
Proinsulin

## Bone Metabolism

Calcitonin  
Osteocalcin  
25-OH Vitamin D  
Intact PTH  
\* $\beta$ -CrossLaps ( $\beta$ -CTX)  
\*total P1NP

\* Available soon

# MAGLUMI® X3

## Fully Automatic Chemiluminescence Immunoassay (CLIA) System

### Technical Specifications

#### Chemiluminescence Principle

- Direct chemiluminescence
- Magnetic microbeads separation, fast and efficient
- Key technology of ABEI-label and magnetic microbead separating enhances stability and sensitivity of MAGLUMI® reagents

#### Main Features

- Throughput: up to 200 tests/hour
- Space occupied <math>0.68\text{m}^2</math>, throughput of per unit area >math>294\text{T/h/m}^2</math>
- 24 hours ready to use

#### Reagent Loading

- 20 reagent positions on board with no-pause loading/unloading
- 24-hour refrigerated, constant at  $10^{\circ}\text{C}$  ( $\pm 0.3^{\circ}\text{C}$ ) in reagent area
- RFID reading all information of reagent
- Indicator light available, real-time status monitoring

#### Sample Loading

- Up to 72 sample positions
- No-pause loading/unloading, STAT available
- Indicator light available, real-time status monitoring
- Barcode recognition: Barcode reader available
- Types of sample tube: micro sample cup, blood collection tube, plastic test tube
- Sample editing mode: LIS connection, barcode recognition, manual editing

#### Reagent Features

- MAGLUMI® reagent kits are compatible with all MAGLUMI® analyzers
- Calibrators and internal controls included
- 10 points master curve, 2-points re-calibration
- RFID tag storing all informations of reagents
- Storage temperature  $2^{\circ}\text{C}$ - $8^{\circ}\text{C}$

#### Reaction System

- Single reaction cup with five-sided heating at  $37.0^{\circ}\text{C}$  ( $\pm 0.3^{\circ}\text{C}$ ) incubation, ensuring accurate, quick and uniform incubation
- Incubating 80 cuvettes simultaneously
- Non-contact mixing with different modes before incubation, ensuring sufficient reaction

#### Sampling System

- Titanium needle for crushproof
- Teflon coated for carry-over prevention
- Liquid level detection, clot detection, prevent bubble interference, collision detection, wash automatically

#### Consumable System

- Single cup design, no-pause loading/unloading
- Starters: 230mL of Starter 1 and Starter 2
- System liquid: Up to load 10L at one time
- Indicator light monitors real-time consumable quantity

#### Washer

- 4-step washing technology
- Intelligent magnetic field design ensure excellent separation
- Three independent non-contact magnetic microbeads vortex cleaning technology
- Resuspension washing mode efficiently reduces non-specific binding and improves sensitivity
- Independent washing unit for waste liquid needle can completely wash without residue

#### Measuring Chamber

- High-sensitivity and low-noise photomultiplier tube (PMT)
- Heating and cooling bidirectional temperature control function included
- Original light pollution isolation technology prevents light pollution and ensures accurate and reliable results

#### Software System

- Friendly software UI design, humanized instrument operation procedure
- Real-time status display for each test procedure
- Self-recovery function make sure the analyzer running smoothly without unnecessary manual intervention
- Bidirectional Snibelis™ connection by TCP/IP and COM Port
- Supporting SnibeLinker™ remote diagnosis system

#### Dimensions & Weight

- Dimensions: 78 (h) × 90 (w) × 75 (d) cm
- Weight: 127 kg

MAGLUMI®, Snibelis™ and SnibeLinker™ are trademarks of Snibe.  
All other product names and trademarks are the property of their respective owners.



Shenzhen New Industries Biomedical Engineering Co., Ltd.  
(SNIBE Co., Ltd)

No.23, Jinxiu East Road, Pingshan District, 518122 Shenzhen,  
P.R. China

Tel: +86 755 26501514  
Email: sales@snibe.com

Fax: +86 755 26654850  
Web: www.snibe.com

Follow us on



## AC6000 Series



### Full-auto Glycated Hemoglobin HbA1c Analyzer

#### Specifications

Model	AC6601 Full-auto glycated hemoglobin HbA1c Analyzer
<b>Basic Parameters</b>	
Testing Method	Chromatography/ion-exchange liquid chromatography
Testing Item	Glycated hemoglobin HbA1c
Testing Scope	3% - 18%
Testing Parameters	Precision (CV) ≤3%, Accuracy (V) ≤1.50%
Testing Time	4 minutes (including reporting)
Sample type	Venous blood (EDTA anticoagulation), Peripheral finger blood
Sample volume	5 μl whole blood for each sample
<b>Functional Parameters</b>	
Photometer	415nm LED integral flow colorimeter
Sampling mode	Full-auto 25-position sample turntable (20 sample positions, 2 quality control positions, 1 emergency position, 1 cleaning position, 1 zero position)
Testing mode	Auto test for batch, single test for emergency
Calibration	Optional selection of 2 points or 4 points calibration, optional selection of manual or automation, Equipped with 1 set of weight high and low level calibrators
Thermostatic control	Column and reagent thermostatic control: 25 ± 0.5°C
Reagent matching	One piece of column and reagent consumable/300 pieces of testing samples
Reagent allocation	Each standard package includes A, B, C, D eluent, calibrator, hemolytic agent, pump tubing
<b>Output Parameters</b>	
Display	320x240 LCD Graphic Display, displaying real-time testing curve
Printer	58mm Thermal printer, printing testing curve and report
Report Output	IFCC concentration value, NGSP area percentage, ADAG average glucose
Data storage	1000 test report (including testing curve)
Communication interface	RS232 communication interface, connecting to HIS/LIS system
<b>Working Parameters</b>	
Power Supply	AC220V ± 22V 50HZ 215VA
size	408mm x 390mm x 230mm
Weight	9kg
Working Environment	Temperature 12°C-30°C, Relative humidity ≤85%

# AC6601 HbA1c

## Viewable Accuracy

new

•Faster • More accurate • More convenient•

Batch testing No wait Inside heating Precision control One step testing





# Four new comprehensive upgrades

## Faster More accurate More convenient

upgrade

### Inside heating

No pre-heating step outside, operate easily.

New

### Durable pump tubing

Maintenance free micro pump, accurate sample volum, accurate results, imported high performance pump tubing, using life is for 1000 tests. Reduce replacement frequency, easier operation.

upgrade

### Filter for avoiding blockage

With filter device for sample, avoid blockage, reduce failure rate greatly.

New

### Pressure monitor

With pressure testing device, monitor pressure at all times, results are more credible.

## HbA1c Accuracy is the main key

It is the world consensus that glycated hemoglobin (HbA1c) becomes the important diagnostic standard in diabetes. HbA1c > 6.5% becomes the starting point for diabetes, and HbA1c > 6.0% and < 6.5% is the main attribute for diabetes' high risk group. The accuracy of HbA1c is utmost important for diabetes screening, diagnosis, blood-glucose control and therapeutic effect monitoring. This interval of HbA1c is very narrow which is only 0.5%. For this reason, the accuracy of analyzer for glycated hemoglobin is critical.

## Real-time Chromatogram True Reliability

AC6601 adopts the display technology with real-time summarized chromatogram, carries out the basic principle of analysis instruments and appears the real original data. Discard modification of test results by means of software during late stage, it is ensured that the whole process of testing can be redisplayed truly, and provide credible results to clinic.

## AC6601 Accuracy + High Performance | Low Price

AC6601 adopts the principle of test methodology same as high pressure liquid chromatography HPLC and uses low pressure liquid chromatography in combination with the advanced test and analysis technology. It's results fully meet the testing requirements for glycated hemoglobin HbA1c in lab and provide accurate test report for clinical diagnosis. The price is much lower than the price of HPLC and reasonable consumable prices enable AC6601 has high performance/Low price. AC6601 is specially suitable for the needs to test glycated hemoglobin HbA1c in all kinds of hospitals.

## 6 Items of Test Technology Ensuring Accurate Results

### ■ Accurate Principle of Methodology

Adopts classical and accurate principle of methodology-ion exchange liquid chromatography. It is the gold standard of HbA1c analysis, and it is the only analysis method to really separate HbA1c directly by measuring the piecemeal absorbance through continuous test on line, and obtain the correct area percentage with integration.

### ■ Accurate separation 4-gradient elution

The novel 4-gradient elution for HbA1c can separate accurately glycated hemoglobin with 4-gradient elution of corresponding concentration reagent aiming at HbA1c instead of routine elution process produced by high and low concentration mixture.

### ■ High Separation Liquid Chromatographic Column

High separation liquid chromatographic column made of imported resin with volume of  $\Phi$  9mm x 45mm and weight of 2.5g which is 15-20 times greater than general micro chromatographic column. High efficiency chromatographic column for 300 tests ensure the accuracy of test results.

### ■ High Sensitivity 415nm LED Integral Photometer


High sensitivity 415nm LED integral photometer has the characteristics of correct wavelength, stable light source, full aluminum alloy structure, high anti-interference performance, multi-lens focusing, micro cuvettes and high sensitivity. It can record accurately analysis curve.

### ■ Supply with Original Calibrator

The international standard value tracing is adopted. The authoritative reference material for quality control is used to transfer the numeric value. Each set of reagent is supplied with two sets of calibrator for proper calibration so as to ensure that the test results are correct and reliable and thoroughly to avoid individual error caused by factor calibration.

### ■ Precision Chromatographic Column and Thermostatic Apparatus

Precision chromatographic column and reagent thermostatic control apparatus ensure that Chromatographic column and reagent are not affected by environment temperature and effectively guarantees the repeatability and correctness of the test results.

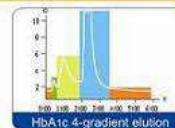

## 6 Items of Human Oriented Design, Convenience, Speed and Reliability

### ■ Real-time summarized chromatogram, Intelligence process detection

The advanced embedded micro processing system + intelligence control software can really display testing curve and the real redisplay of testing process and can monitor and alarm for test results, absorbency, signal potential, peak time and reagent consumption. It can be operated with high proficiency.

### ■ Full-auto 25-position Sample Turntable, Optional Selection of Batch Test or Emergency Test

With full-auto 25-position sample turntable, it is unnecessary to use complicated rank sampling device and built-in hemolysis device. Batch test can be automatically carried out. Emergency test can be carried out at any time.

### ■ Full-open Structure, perfect flow path, low failure and easy maintenance

Full-open structure, Perfect solenoid valve flow path, it is unnecessary to use complicated sample mechanical 6-way and rotating distribution control valve. It is reliable for use and easy for maintenance.

### ■ 5 $\mu$ l whole blood, for both labs and clinics



5 $\mu$ l whole blood, both venous blood and peripheral finger blood can be tested, suitable for labs for batch tests and also suitable for clinics for providing test results immediately.

### ■ With gas solution and bubble-removal technology, the error can be eliminated overall

Since the instrument is equipped with reagent solvability gas eliminating device and applies the cuvette bubble auto detection eliminating technology, the air bubbles which affect the test results can be eliminated overall without the complicated degassing device.

### ■ Saccharification concentration, area percentage and average glucose can be synchronously reported

The instrument can output the data including IFCC concentration value, NGSP area percentage and ADAG average glucose on the test report synchronously to meet the requirement of world standardization. It may memorize 1000 test curve reports and is equipped with RS232 communication interface and can connect directly to HIS/LIS system.

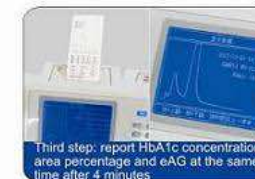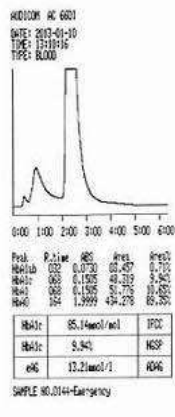

**Tumor  
Markers**

**Hormone  
Markers**

**Diabetes  
Markers**

**Cardiac  
Markers**

**Infection  
Markers**

# fluorecare<sup>®</sup>

Immunofluorescence Quantitative Platform

# fluorecare®



Microprofit fluorecare® Immunofluorescence Quantitative Platform uses immunochromatographic technology, for the diagnosis of common human diseases. Based on specific antigen-antibody reaction, the triggered fluorescence from the fluorochrome-labeled complexes at the test line is collected and automatically calculated. For certain analyte, including serum, plasma, urine or other sample, fluorecare® is able to measure the concentration of various biological markers including tumor markers, hormone markers, infection markers, cardiac markers and diabetes markers with excellent accuracy and easy operation.


## ◆ Test items

- ▶ Tumor  
CEA, AFP, PSA, f-PSA, Cyfra21-1, CA125, CA19-9, PG I/II, Ferritin, Transferrin
- ▶ Hormone  
β-HCG, LH, FSH, PRL, Testosterone, Progesterone, AMH, TSH, T3, T4, 25-OH Vitamin D, Cortisol
- ▶ Diabetes  
HbA1c, Insulin, C-Peptide, MAU, CYS-C
- ▶ Infection  
PCT, CRP, SAA
- ▶ Cardiac  
cTnI, MYO, CK-MB, NT-proBNP, D-Dimer

## ◆ Advantages

- POCT (Point-of-care testing)
- Room temperature shipping and storage
- Full platform, multiple parameters
- Inner and independent quality control

## ◆ Technique Principle

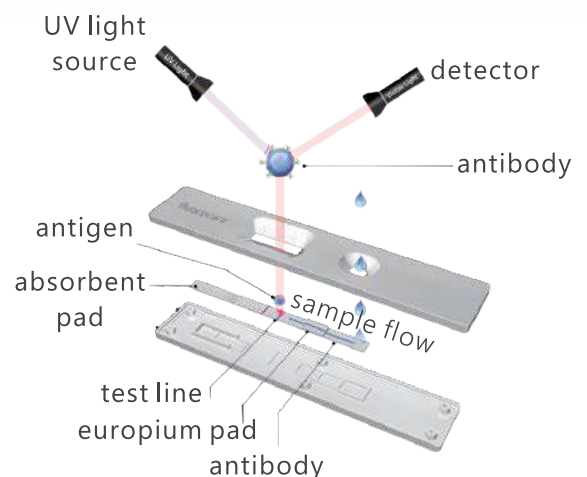

- One-step rapid test



# fluorecare<sup>®</sup> MF-T1000

## Immunofluorescence Analyzer

- ▶ **No specialized facilities or personnel needed.**  
Unlike traditional workstation, fluorecare POCT platform does not require specialized facilities. Easy to learn, easy to use.
- ▶ **Automatic calibration and comprehensive parameters.**  
In-built automatic calibration saves time and maintains accuracy.
- ▶ **Intended use in community hospitals and primary clinics.**  
Include primary clinics, small laboratories, community hospitals, medical facilities in rural areas and emergency department in hospitals.
- ▶ **Abundant connectivity options.**  
4 USB ports, 1 LAN port, 1 LIS port and 1 PC connection port.


# fluorecare<sup>®</sup> MF-C100

## Temperature Control Instrument





- ▶ **Precise temperature control.**
- ▶ **Automatic timing function.**
- ▶ **A perfect match for fluorecare<sup>®</sup> MF-T1000 and test cassettes.**



## ◆ Workflow

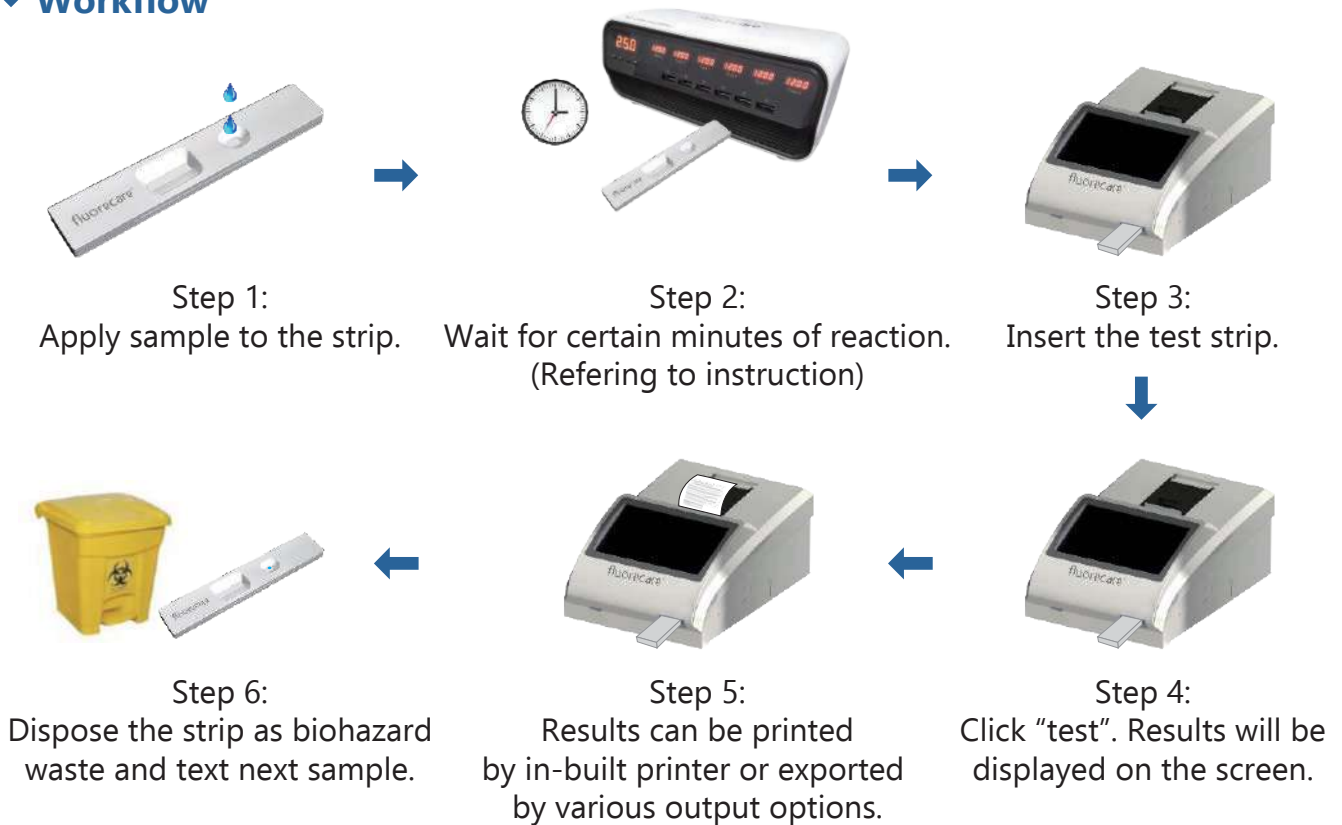

## ◆ High consistence of fluorecare® Immunofluorescence Quantitative Platform with leading product manufacturer.

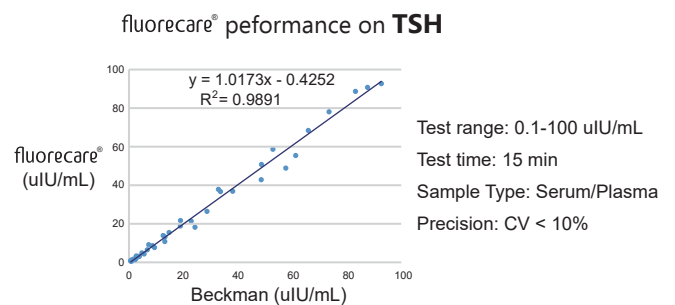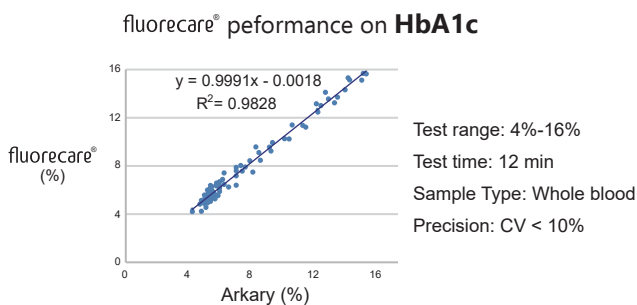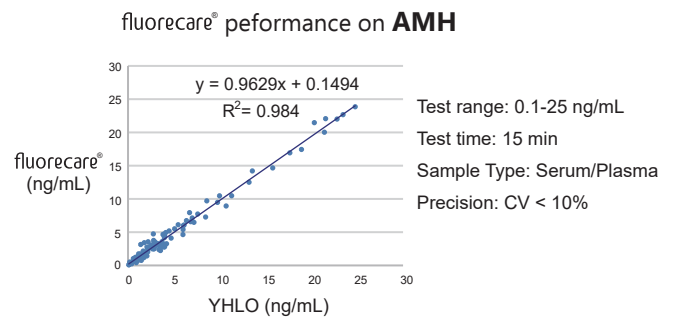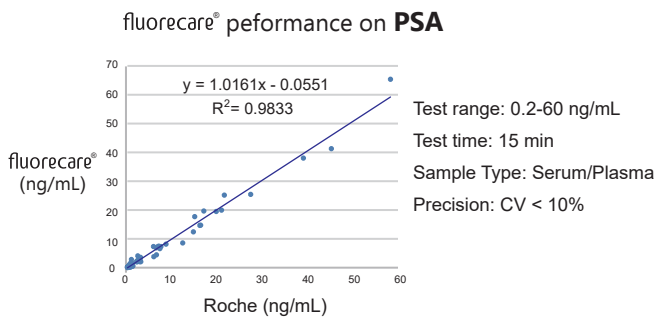

Product Class	Markers	Specimen	Reaction Time	Measuring Range	Qualification	Ongoing
<b>Tumor</b>	CEA	S/P	15 min	1-200 ng/ml	CE	
	AFP	S/P	15 min	2-300 ng/ml	CE	
	PSA	S/P	15 min	0.2-60 ng/ml		
	f-PSA	S/P	15 min	0.1-25 ng/ml		
	Cyfra21-1	S/P	15 min	1-200 ng/ml	CE	
	CA125	S/P	15 min	1-500 U/ml	CE	
	CA19-9	S/P	15 min	5-300 U/ml	CE	
	CA15-3	S/P	15 min	1-300 U/ml	CE	▲
	PG I/II	S/P	10 min	PG I: 2-800 ng/ml/PG II: 1-100 ng/ml	CE	
	Ferritin	S/P	15 min	10-1000 ng/ml	CE	
	Transferrin	S/P	15 min	0.25-5 g/L	CE	
	Pro-GRP	S/P	15 min	4-500 pg/ml	CE	▲
	HE-4	S/P	15 min	50-1500 pmol/ml	CE	▲
	NSE	S/P	15 min	1-200 ng/ml	CE	▲
	SCCA	S/P	15 min	0.3-100 ng/ml	CE	▲
<b>Hormone</b>	β-HCG	S/P	15 min	5-200000 mIU/ml	CE	
	LH	S/P	15 min	1-100 mIU/ml	CE	
	FSH	S/P	15 min	1-100 mIU/ml	CE	
	PRL	S/P	15 min	1-100 ng/ml	CE	
	Testosterone	S/P	15 min	0.1-15 ng/ml	CE	
	Progesterone	S/P	15 min	0.8-40 ng/ml	CE	
	AMH	S/P	15 min	0.1-25 ng/ml	CE	
	TSH	S/P	15 min	0.1-100 uIU/ml	CE	
	T3	S/P	10 min	0.5-12 nmol/L	CE	
	T4	S/P	10 min	10-300 nmol/L	CE	
	ft3	S/P	15 min	1-50 pmol/ml	CE	▲
	ft4	S/P	10 min	3-77 pmol/L	CE	▲
	25-OH Vitamin D	S/P	10 min	3-100 ng/ml	CE	
	Cortisol	S/P	15 min	0.5-75 µg/dL	CE	
<b>Diabetes</b>	HbA1c	WB	12 min	4%-16%	CE	
	Insulin	S/P	15 min	2-300 mIU/L	CE	
	C-Peptide	S/P	15 min	0.1-30 ng/ml	CE	
	MAU	Urine	15 min	5-200 mg/L	CE	
	CYS-C	S/P	3 min	0.1-8 mg/L	CE	
<b>Infection</b>	PCT	S/P	15 min	0.1-50 ng/ml	CE	
	CRP	S/P/WB	3 min	0.5-150 mg/L	CE	
	SAA	S/P	5 min	5-150 mg/L	CE	
<b>Cardiac</b>	cTnI	S/P	15 min	0.1-50 ng/ml	CE	
	MYO	S/P	15 min	1-500 ng/ml	CE	
	CK-MB	S/P	15 min	1-80 ng/ml	CE	
	NT-proBNP	S/P	15 min	100-30000 pg/ml	CE	
	D-Dimer	S/P	10 min	0.1-10 mg/L	CE	

▲ Under development

Rev.: 2019/07/01

fluorecare<sup>®</sup>, creating a **NEW** impact on point-of-care testing.



Microprofit Biotech Authorized Distributor

Distributed by:

**ImTec**

**DIAGNOSTICS NV**

Contact: [info@imtec.be](mailto:info@imtec.be) | +32(0)3-2251139



# URIT

High reliability & automatization with convenience

## URIT-500B Urine Analyzer



Global Diagnostics Supplier



# URIT-500B Urine Analyzer

FDA 510 ( K ) No.K082811 Approved



## Specifications

**Parameters:** Leukocyte, Ketone, Nitrite, Urobilinogen, Bilirubin, Protein, Glucose, Specific Gravity, Blood, pH, Ascorbic Acid, Microalbumin, Creatinine, Calcium

**Measurement Principle:** Dual wavelength reflectance photometry

**Reaction time:** 65 seconds

**Throughput:** 520 samples /hour

**Display:** LCD touch screen

**Memory Capacity:** 9999 results

**Strip:** URIT 10G/11G/14G urine reagent strip.

**Interface:** Serial port(RS232), parallel port (For external printer), PS/2(For barcode reader),USB Port(For data out put)

**Work condition:** Environmental Temperature: 15-30°C,  
RH: 20%~80%,Lumination:=2000 lx

**Power supply:** 100~240V AC,50/60Hz

**Dimension:** 390mm x 340mm x 290mm(Depth×Width×Height)

**Weight:** 6.5kg

## Features & Functions

- ▣ Friendly operation with touch screen
- ▣ Easy maintenance design with auto waste collection
- ▣ High test speed, ideal for extensive sample laboratory
- ▣ Fully automatic operation for ready-to-test
- ▣ Internal calibration for sensitivity of strips
- ▣ Rs232 serial interface to host PC/LIS, Barcode reader optional
- ▣ URIT UQ-10,UQ-11 and UQ-14 Urinalysis Control



## URIT 10G、11G、14G urine reagent strip

14G with Micro albumin, Creatinine and Calcium.

## URIT Medical Electronic Co., Ltd.

Head Office

Add.: No. D-07 Information Industry District,

High-tech Zone, Guilin, Guangxi 541004,

P. R. China

Tel: +86-773-2288586

Fax: +86-773-2288560

Shenzhen Office

Add.: Room E.J, 24th Floor, New Baohui Building, Nanhai Road,

Nanshan District, Shenzhen, P.R. China

Tel: +86-755-26050853

Fax: +86-755-86036150

E-mail: export@uritest.com service@uritest.com [Http://www.urit.com](http://www.urit.com)



# GH-900Plus HbA1c Analyzer

## TECHNICAL SPECS

Methodology	High - Performance Liquid Chromatography (HPLC)
Test Modes	Variant Mode (A1c\E\D\S\C)
Test Range	3% - 18%
Precision	CV ≤ 1.5%
Test Speed	130 Secs / Test
Sample Type	Venous Blood, Finger Peripheral Blood, Lyophilized Whole Blood 10μL(whole blood), 400μL (Diluted blood)
Auto Sample Station	5 Positions/Rack
Photometer	415nm + 500nm Detector
Filter	≥ 400T
Display	10.1 " TFT True Color LCD Touch Screen
Software	Linux Software with Self - Diagnosis to Monitor and Detect System Errors
Reagent Kit	Eluent A, Eluent B, Hemolysin L, Calibrator; QC Material
Information Input	Scanner or Touch Keypad
Storage	4000 Sample Results
Connection	USB, LAN, LIS Compatible
Printer	Thermal Printer
Operation	Temperature 10 ~ 30 °C ( 50 ~ 86 °F )
Humidity	≤ 80%
Power	AC 100-240V 50/60HZ 120VA
Dimensions	450mm(L)*360mm(W)*540mm(H)
Weight	32.8KG



# Lifotronic

— Caring for Better Life —



## GH-900Plus HbA1c Analyzer

### Lifotronic Technology Co., Ltd.

**Address:** Unit A,4th Floor, Building 15, Yijing Estate, No. 1008 Songbai Road, Nanshan District, Shenzhen City, Guangdong Province, 518055, P.R.China

**Tel:** 86-755-29060026 **Fax:** 86-755-29060036

**Email:**inter-marketing@lifotronic.com **Web:** en.lifotronic.com

[en.lifotronic.com](http://en.lifotronic.com)

# GH-900Plus HbA1c Analyzer



Lifotronic fully automated GH-900Plus HbA1c Analyzer supports fast HbA1c results output in 130 seconds without Hb variant interference. It provides the outstanding solution for reliable diabetic monitoring. No sample preparation and walk-away operation.

### HPLC Technology – Gold Standard Methodology

- NGSP and IFCC Certified
- HbA1c Results within 130 Seconds
- 5 samples loading capacity, suitable for medium to small labs

### Fully Automated - To Minimize Operation Steps

- No sample preparation
- Fully automated system cleaning after test
- External barcode scanner for sample identification

### Precise and Reliable – To Serve You Consistently

- Inter measuring  $CV \leq 1.5\%$  & Intra measuring  $CV's \leq 3\%$
- Superior quality chromatographic resolution to eliminate interferences

### Dual Wavelength Detection – To Avoid Interference



- To avoid the reagent peak interference
- More anti-interference abilities, the mutation factor interference to the peak can be easily counteracted
- To eliminate the nonspecific absorption of hemoglobin

### Degasser – For Better Result Accuracy

- More Stable Pressure, More Accurate Flow Rate
- To Reduce Background Absorption and Improve Detection Sensitivity
- To Improve the Separation Effect of Column and Prolong Its Lifetime


### Compact Size – To Minimize Space Requirements

- Small Footprint Reduces Bench Space Needed
- Most compact fully automated HPLC system

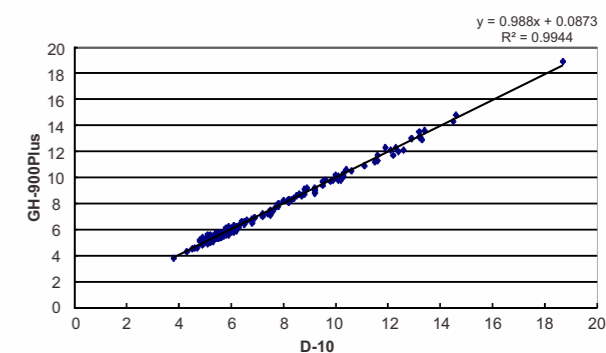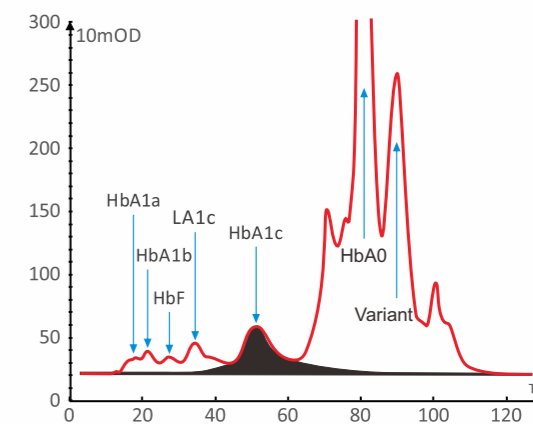

Correlation between GH-900Plus and D-10

### The Elements of Hemoglobin

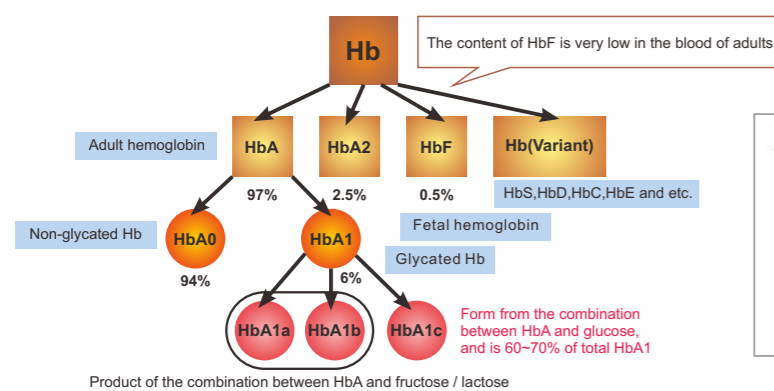

IFCC



NGSP

### Gold Standard of Diabetes Diagnose

Glycosylated hemoglobin (HbA1c) is widely recognized as a Gold Standard to monitor diabetes, which can indicate the average plasma glucose concentration over 8 ~ 12 weeks.

### HPLC Methodology

High-Pressure Liquid Chromatography (HPLC), to separate HbA1c directly with measuring the absorbance points continually to form chromatogram. Using normal distribution curve fitting auto-iterative algorithm to get precise HbA1c testing result, excluding interference of variant and unstable hemoglobin like LA1c. Standard Analysis Mode will report HbA1a, HbA1b, HbF, LA1c, HbA1c, HbA0 peak areas and percentage. And the result also includes IFCC, NGSP and ADAG value for diverse client needs.

# H8 Hemoglobin Analyzer (HPLC)

## TECHNICAL SPECS

<b>Methodology</b>	High-Performance Liquid Chromatography (HPLC)
<b>Test Modes</b>	Fast Mode, Variant Mode
<b>Test Range</b>	3% - 18%
<b>Precision</b>	CV ≤1.5%
<b>First Sample Result</b>	7.5 Mins
<b>Test Speed</b>	1.5Mins/sample for Fast Mode, 2.2 Mins/sample for Variant Mode
<b>Sample Type</b>	Venous Blood, Finger Peripheral Blood, Lyophilized Whole Blood Peripheral Blood, 500ul(150 Dilution Ratio)
<b>Auto Sample Station</b>	10 Positions
<b>Photometer</b>	415nm+500nm LED, 20000 Hours Life Span
<b>Chromatography Column</b>	Available Tests ≥800T
<b>Filter</b>	≥400T
<b>Display</b>	10.1"TFT True Color LCD Touch Screen
<b>Software</b>	Embedded System with Self-Diagnosis to Monitor and Detect System Errors
<b>Reagent Kit</b>	Eluent A, Eluent B, Eluent C, Hemolysin, Calibrator, QC Material (Weight Sensor±1%)
<b>Information Input</b>	Scanner or Touch Keypad
<b>Storage</b>	4000 Sample Results
<b>Connection</b>	USB, LAN, LIS Compatible
<b>Printer</b>	Thermal Printer and External Laser Printer
<b>Operation</b>	Temperature 10~30°C (~°F)
<b>Humidity</b>	≤85%
<b>Power</b>	AC 100-240V 50/60HZ 120VA
<b>Dimensions</b>	580mm×500mm×520mm (22.8"H×19.7"W×20.5"D)
<b>Weight</b>	50kg (110lbs)
<b>Barcode Scanner</b>	QC Curve



**Lifotronic Technology Co., Ltd.**

Address: 4th Floor, Building 15, 1008 Songbai Road, Nanshan District, Shenzhen, 518055, China

Tel: 86-755-29060197 Fax: 86-755-29060036 Email: Inter-sales@lifotronic.com

Web: en.lifotronic.com



# H8 Hemoglobin Analyzer (HPLC)

[en.lifotronic.com](http://en.lifotronic.com)

# H8 Hemoglobin Analyzer (HPLC)



The Fully automated Lifotronic H8 Hemoglobin Analyzer offers the fast throughput of HbA1c results in 130 seconds, with Hb variant detection, providing the outstanding solution for quick and reliable diabetic monitoring. No sample preparation and very little hands-on time by the operator is required for the H8 Analyzer.




## Gold Standard of Diabetes Diagnose


Glycosylated hemoglobin (HbA1c) is widely recognized as a Gold Standard to monitor diabetes, which can indicate the average plasma glucose concentration over 8~12 weeks.

## HPLC Methodology

High-Pressure Liquid Chromatography (HPLC), to separate HbA1c directly with measuring the absorbance points continually to form chromatogram. Using normal distribution curve fitting auto-iterative algorithm to get precise HbA1c testing result, excluding interference of variant and unstable hemoglobin like HbF. Standard Analysis Mode will report HbA1a, HbA1b, HbF, HbA1c, LA1c, HbA0 peak areas and ratio. And the result also includes IFCC, NGSP and ADAG value for diverse client needs.

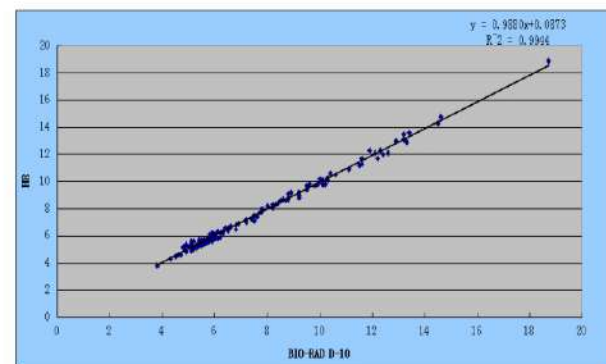


Correlation between H8 and a famous HPLC method

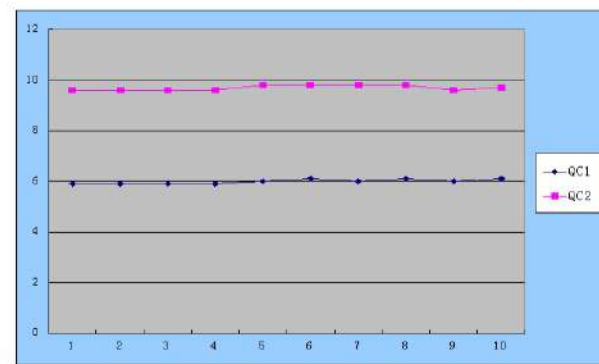

Precision Study

## The Elements of Hemoglobin

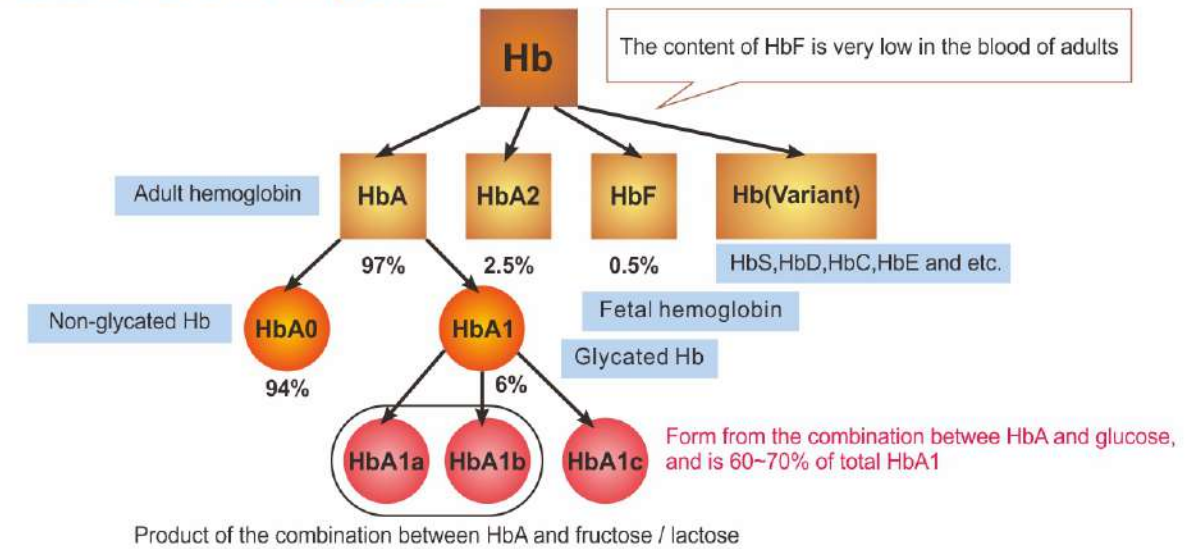

## HPLC Technology – Gold Standard Methodology

- NGSP and IFCC certified
- HbA1c results in 130 seconds
- Fully Automated - To Minimize Operation Hassles
  - Primary Tube Sampling with Cap Piercing
  - Fully Automated Start-up, Maintenance and Shutdown
  - Barcode Scanner for Sample identification
- Precise and Reliable – To Serve You Consistently
  - HbA1c Inter & Intra Measuring CV's  $\leq 3\%$  to Enable Exceptional Result Management
  - Superior Quality Chromatographic Resolution to Eliminate Interferences

## Dual Wavelength Effect

- To Avoid the Reagent Peak Interference
- More Anti-interference Capabilities, the Interference of the Mutation Factor to the Peak Can Easily be Counteracted
- To Eliminate the Nonspecific Absorption of Hemoglobin

## Dual Wavelength Effect

- More Stable Pressure, More Accurate Flow Rate
- To Reduce Background Absorption and Improve Detection Sensitivity
- To Add the Separation Effect of Column and Prolong Its Lifetime

## Compact Size - To Minimize Space Requirements

- Small Footprint Reduces Bench Space Needed

